D13

Europäisches Patentante
Europaan Patent Office
Office européen des brevets

19 Publication number: 0 2

0 268 237

ூ

EUROPEAN PATENT APPLICATION

2) Application number: 87116861.3

② Date of filing: 16.11.87

(9) Int. CI.4: G01N 1/10 , G01N 35/00 , G01F 11/02

Driority: 17.11.86 US 931476

Date of publication of application: 25.05.88 Bulletin 88/21

Designated Contracting States: AT BE CH DE ES FR GB GR IT LI LU NL SE 2 Applicant: ABBOTT LABORATORIES

Abbott Park Illinois 60064(US)


nventor: Hayes, Donald J.
2012 Tampicko Drive
Plano Texas 75075(US)
inventor: Wallace, David B.
9929 Wood Forest
Dallas Texas 75243(US)
inventor: Verlee, Donald J.
553 Drake Street
Libertyville Illinois 60048(US)
inventor: Houseman, Kenneth R.
1520 S. Main Street
Racine Wisconsin 53403(US)

Representative: Modiano, Guido et al MODIANO, JOSIF, PISANTY & STAUB Modiano & Associati Via Meravigii, 16 i-20123 Milan(IT)


Apparatus and process for reagent fluid dispensing and printing.

A system for printing and dispensing chemical reagents in precisely controlled volumes onto a medium at a precisely controlled location. A jetting tube, comprising an orifice at one end and a fluid receiving aperture at the other end, is concentrically mounted within a cylindrical piezo-electric transducer. The fluid receiving aperture is connected to a reservoir containing a selected reagent by means of a filter. The reservoir is pressurzed by a regulated air supply. An electrical signal of short duration is applied to the transducer. The pulse causes the transducer and the volume defined by the jetting tube to expand, thereby drawing in a small quantity of reagent fluid. The cessation of the pulse causes the transducer and the volume of the jetting tube to de-expand, thereby causing at least a substantially uniformly sized droplet of reagent fluid to be propelled through the orifice. The droplet may be directed to impact a printing medium or collected in a dispensing recepticle.



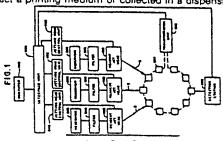

Xerox Copy Centre

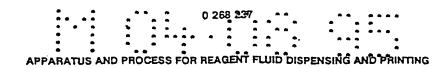

BACKGROUND OF THE INVENTION

The present invention relates to an apparatus and process for dispensing and printing reagent fluids, wherein a transducer is used to propel small quantities of the fluid towards a positioned target.

Diagnostic assays often require systems for metering, dispensing and printing reagent fluids. In the case of metering and dispensing, such systems comprise both manual and automatic means. For purposes of practicality, the present background discussion will focus on the methods of metering and dispensing 100 micro-liter volumes or less.

The manual systems of metering and dispensing include the glass capillary pipet, the micro-pipet, the precision syringe; and weighing instruments. The glass capillary pipet is formed from a precision bore glass capillary tube. The pipet typically comprises a fire blown bulb and a tubular portion fire drawn to a fine point. Fluid is precisely metered by aspirating liquid through the tube into the bulb to a predetermined level indicated by an etched mark. The fluid may then be dispensed by blowing air through the tube.

The micro-pipet typically comprises a cylinder and a spring loaded piston. The travel of the piston is precisely determined by a threaded stop. The distance the piston travels within the cylinder and the diameter of the cylinder define a precise volume. The fluid is aspirated into and dispensed from the micro-pipet in precise quantities by movement of the piston within the cylinder.

The precision syrings generally comprises a precisely manufactured plunger and cylinder with accurately positioned metering marks. The fluid is introduced into and dispensed from the syringe by movement of the plunger between the marks.

Weighing techniques for dispensing fluids often simply involve weighing a quantity of fluid. The density of the fluid may then be used to determine the fluid volume.

Exemplary automatic metering and dispensing systems include the precision syringe pump; the peristaltic pump; and the high performance liquid chromatography (HPLC) metering valve. The precision syringe pump generally comprises a precision ground piston located within a precision bore cylinder. The piston is moved within the cylinder in precise increments by a stepping motor.

The peristaltic pump comprises an elastomeric tube which is sequentially pinched by a series of rollers. Often the tube is placed inside a semi-circular channel and the rollers mounted on the outer edge of a disc driven by a stepping motor. The movement of the rollers against the tubing produces peristaltic movement of the fluid.

The HPLC metering valve comprises a defined length of precision inner diameter tubing. The fluid is introduced into the define volume of the tubing with the valve in a first position and then dispensed from the tubing when the valve is placed in a second position.

All of the above metering and dispensing systems have the disadvantage that the volumes dispensed are relatively large. Furthermore, these systems are also relatively slow, inefficient and comprise precision fitted components which are particularly susceptible to wear.

The printing of reagent fluids is frequently required in the manufacture of chemical assay test strips. Selected reagents are printed in a desired configuration on strips of filter paper. The strips may then be used as a disposable diagnostic tool to determine the presence or absence of a variety of chemical components.

'Generally, to perform a chemical assay with a test strip, the strip is exposed to a fluid or a series of fluids to be tested, such as blood, serum or urine. In some instances, the strip is ninsed and processed with additional reagents prior to being interpreted. The precise interpretation depends on the type of chemical reactions involved, but it may be as simple as visually inspecting the test strip for a particular color change.

The manufacture of test strips generally involves either a manufacturing process or a blotting process. The blotting process is the simplest manufacturing method and permits most reagents to be applied without modification. A disadvantage of this process is that it is difficult to blot the fluids onto the test strip with precision.

The printing process will often involve any of three well known methods: silk screening: gravure: and transfer printing. The silk screening of reagents generally involves producing a screen by photographic methods in the desired configuration for each reagent to be printed. The screen is exposed under light to a preselected pattern and then developed. The areas of the screen which are not exposed to light, when devel oped, become porous. However, the areas of the screen which have been exposed to light remain relatively nonporous. The screen is then secured in a frame and the test strip placed below. The desired

reagent fluid, Specially prepared to have a high viscosity, is spread over the top side of the screen. The reagent passes through the porous areas of the screen and onto the test strip. The test strip is then subjected to a drying process, specific to each reagent. Once the test strip is dry, it may be printed again using a different screen, pattern and reagent.

The gravure method of printing reagents comprises coating a metal surface with a light sensitive polymer. The polymer is exposed to light in the desired predetermined pattern. When developed, the polymer creates hydrophilic and hydrophobic regions. The reagent is specially prepared such that when applied to the metal it will adhere only to the hydrophilic regions. After the specially prepared reagent is applied, the test strip is pressed against the metal and the reagent is transferred from the metal to the test strip.

The transfer printing method comprises transferring the reagents from a die to the test strip in the desired pattern. The die is made with the appropriate pattern on its surface and then coated with the desired, specially prepared reagent. A rubber stamp mechanism is pressed against the die to transfer the reagent in the desired pattern from the die to the rubber stamp. The rubber stamp is then pressed against the test strip to transfer the reagent, in the same pattern, to the test strip.

Each of the above-mentioned reagent printing techniques has significant disadvantages. The most common disadvantage is the requirement that the reagents must be specially prepared. Additionally, if a variety of reagents are to be printed onto a single test strip, the strip must be carefully aligned prior to each printing. This alignment procedure increases the cost and decreases the throughput of the printing process. Moreover, a special die or screen must be produced for each pattern to be printed. A further disadvantage arises in that the above printing methods are unable to place reproduceable minute quantities of reagent on the test strip.

It is an object of the present invention to provide a printing and dispensing method and apparatus which avoids these disadvantages.

SUMMARY OF THE PRESENT INVENTION

The present invention is directed to a reagent dispensing and printing apparatus and method, wherein the apparatus comprises a transducer operative to eject a substantially uniform quantity of reagent in a precise predetermined direction.


According to one preferred embodiment of the present invention used in dispensing reagent fluids, a jetting tube is concentrically located with a piezoelectric transducer. The jetting tube comprises an orifice at one end and a reagent receiving aperture at the other end. The receiving end of the jetting tube is connected to a filter which is in turn connected to a reservoir containing a selected reagent. A jetting control unit supplies an electrical pulse of short duration to the transducer in response to a command issued by a computer. The electrical pulse causes the volume defined by the jetting tube to expand by an amount sufficient to intake a small quantity of reagent fluid from the reservoir. At the end of the pulse duration, the transducer de-expands propelling a small quantity of the reagent fluid through the onfice and into a fluid recepticle. If desired, additional droplets may be deposited in the recepticle or the recepticle aligned with an additional jetting tube for receiving an additional reagent fluid.

An additional preferred embodiment of the present invention may be used for printing reagent fluids onto a print medium. In this embodiment, the jetting tube is aligned with the printing medium such that the propelled droplet impacts a precise position on the medium. The jetting tube or print medium may then be repositioned and another droplet expelled from the jetting tube. The process may be repeated until a desired configuration of the reagent fluid is printed on the medium.

One advantage of the present invention is that precise minute quantities of reagent fluid may be dispensed or printed in a reproducible manner. Additionally, the method and apparatus may be used to emit droplets of fluids having a wide range of reagent fluid viscosities and surface tensions. The reagents do not in general have to be specially adapted for use with the present invention.

The invention itself, together with further objects and attendant advantages, will best be understood by reference to the following detailed description, taken in conjunction with the accompanying drawings.

55

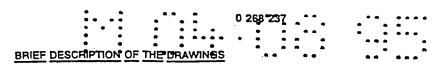

25

35

40

45

FIGURE 1 is a schematic representation of a first preferred embodiment of the present invention showing the use of multiple jetting heads to meter and dispense reagent fluid.

FIGURE 2a is a perspective view of a first preferred embodiment of the jetting head of the present invention.

FIGURE 2b is a cut-away perspective view of the preferred embodiment of Fig. 2a taken along lines 2b-2b with the contact pins removed.

FIGURE 2c is a sectional representation of the preferred embodiment of Fig. 2a taken along lines 2c-2c.

FIGURE 2d is a sectional representation of the preferred embodiment of Fig. 2c taken along lines 2d-2d.

FIGURE 2e is a sectional representation of the jetting tube and transducer of the preferred embodiment of Fig. 2b taken along lines 2e-2e.

FIGURE 3 is a schematic representation of a second preferred embodiment operating in the drop on demand mode as a reagent printing system.

FIGURE 4 is a schematic representation of a third preferred embodiment operating in the continuous mode as a reagent printing system.

FIGURE 5a is a schematic representation of a portion of the jetting head control unit showing the LED strobe circuit.

FIGURE 5b is a schematic representation of a portion of the jetting head control unit showing the high voltage power supply circuit.

FIGURE 5c is a schematic representation of a portion of the jetting head control unit showing the print control circuit.

FIGURE 5d is a schematic representation of a portion of the jetting head control unit showing a portion of the print pulse generator.

FIGURE 5e is a schematic representation of a portion of the jetting head control unit showing an additional portion of the pulse generator.

FIGURE 6a is a perspective view of a second preferred embodiment of the jetting head of the present invention.

FIGURE 6b is an exploded view of the preferred embodiment of Fig. 6a.

FIGURE 7 is a sectional representation of a third preferred embodiment of the jetting head of the present invention.

FIGURE 8 is a sectional view of a symmetrical portion of a fourth preferred embodiment of the jetting head of the present invention.

FIGURE 9 is a graph of the drop mass of the emitted droplets as a function of emission frequency for several fluid viscosities.

FIGURE 10 is a graph of the velocity of the emitted droplets as a function of frequency for several fluid viscosities.

FIGURE 11 is a graph of the total weight of fluid emitted as a function of the number of emitted droplets for a given fluid.

DETAILED DESCRIPTION OF THE PRESENTLY PREFERRED EMBODIMENTS

Turning now to the drawings. Fig. 1 shows a schematic representation of a first preferred embodiment of a reagent dispensing system generally represented as reference numeral 30. The dispensing system 30 comprises a plurality of reagent fluid reservoirs 200, a plurality of filters 300, a plurality of reagent jetting heads 400, a plurality of jetting head control units 500, an interface unit 600, a computer 700, transportation unit 902, a plurality of fluid mixing cells 904 and a detection station 906.

The reservoir 200 holds a selected quantity of reagent fluid for dispensing. The reservoir 200 is maintained at atmospheric pressure by suitable means such as an atmospheric vent. The reagent fluid is transferred from the reservoir 200 through the filter 300 to the reagent jetting head 400. The filter 300 is placed between the reservoir 200 and the jetting head 400 to ensure that any particular foreign matter in the reagent fluid is trapped before entering the jetting head 400.

The plurality of jetting heads 400 and the detection station 906 define a processing path. Each jetting head 400, which is described in detail below, ejects uniformly sized droplets 2 of reagent fluid. The droplets 2 are propelled, with controlled velocity and direction, towards a selecting mixing cell 904 positioned along

0 268,237

the processing path by the traffsportation unit 902. The mixing cells 904 are comprised of non-reactive material and function as minute holding tanks for the dispensed reagent fluid.

The plurality of jetting heads 400, shown in Fig. 1, are positioned sequentially along the processing path. Alternately, some or all of the plurality of jetting heads 400 may be positioned with respect to the transportation unit 902 such that the heads 400 direct the droplets 2 into a selected mixing cell 902 simultaneously.

The jetting heads 400 and the transportation unit 902 are controlled by the computer 700. The computer 700 issues commands to an interface unit 600 which is electrically connected to the transportation unit 902 and to the jetting head control unit 500. The interface unit 600 is of conventional design and is used to control the transfer of information between the computer 700 and the jetting control unit 500. The interface unit 600 is also used to control the transfer of information between the computer 700 and the transportation unit 902.

A first embodiment of the reagent jetting head is shown in Figs. 2a - 2e and generally represented by numeral 400. The jetting head 400 comprises a two piece symmetrical housing 402, 404. The housing 402, 404, when assembled, is adapted to form an orifice aperture 406, an air vent and reagent supply channel 410 and a transducer chamber 403, shown in Fig. 4b. Four screws 408, adapted to respective housing screw apertures 416, hold the housing 402, 404 in an assembled configuration.

The jetting head 400 further comprises a jetting tube 432, a piezo-electric transducer 434 and a reagent fluid supply tube 430. The jetting tube 432 defines a tapered orifice 433 at one end and a fluid receiving aperture 431 at the other end for expelling and receiving fluid, respectively. The piezo-electric transducer 434 is cylindrically shaped and secured concentrically about the mid-region of the jetting tube 432 with epoxy or other suitable means,

The piezo-electric transducer 434, shown in Fig. 2e, defines a first and second end and comprises a section of cylindrically shaped piezeo-electric material 435. An inner nickel electrode 437 covers the inner surface of the cylinder 435. The electrode 437 wraps around the first end of the cylinder 435 a sufficient distance to enable electrical connection external to the cylinder 435.

A second nickel electrode 436 covers the majority of the outer surface of the cylinder 435. The second electrode is electrically isolated from the first electrode 437 by an air gap at the face of the second end of the cylinder 435 and by an air gap on the outer surface of the cylinder 435 near the first end. When an electrical pulse is applied to the first and second electrodes 437, 436 a voltage potential is developed radially across the transducer material 435. The voltage potential causes the radial dimensions of the transducer 435 to change, which causes the volume defined by the transducer 434 to also change.

The jetting tube 432 is positioned in the transducer chamber 403 such that the receiving end 431 extends beyond the rearward end of the transducer 434. The receiving end 431 of the jetting tube 432 is inserted into one end of a reagent supply tube 430. The supply tube 430 is sealingly held to the jetting tube 432 by concentric teeth 412 formed by the housing sections 402, 404. The teeth 412 not only seal the supply tube 430 to the jetting tube 432, but, also, seal the supply tube 430 to the housing 402, 404.

The second end of the supply type 430 passes through the channel 410 and into a reagent reservoir 200. The reservoir 200 contains the reagent fluid to be dispensed by the jetting head 400. As the reagent fluid is dispensed, air is supplied to the reservoir 200 through the channel 410 to prevent the creation of a vacuum in the reservoir 200. The reservoir 200 is releasably attached to the housing 402, 404 and held in place by frictional forces. A reservoir cap 202 is flexibly attached to the reservoir 200 and adapted such that the cap 202 may be used to secure the opening in the reservoir 200 when the reservoir 200 is disengaged from the housing 402, 404.

The position of the jetting tube 432 defines the horizontal plane of the jetting head 400. The jetting tube 432 and the transducer 434 are held in a pre-defined vertical relationship with respect to the housing 402. 404 by means of two upper vertical alignment pins 418 and two lower vertical alignment pins 418. The two upper vertical alignment pins 418 extend horizontally from the housing section 402 into the transducer chamber 403. Similarly, the two lower vertical alignment pins 418 extend horizontally from the housing section 404 into the transducer chamber 403. Each vertical alignment pin 418 is formed integrally with the respective housing sections 402, 404.

The jetting tube 432 and the transducer 434 are held in a predefined horizontal relationship with respect to the housing 402, 404 by means of four horizontal alignment pins 424. Two of the horizontal alignment pins 424 extend horizontally from the housing section 402 approximately midway into the transducer chamber 403. Similarly, two of the horizontal alignment pins 424 extend horizontally from the housing section 404 approximately midway into the transducing chamber 403. Each horizontal alignment pin 424 is formed integrally with the respective housing section 402, 404. The alignment pins 418, 424, sealing teeth 412 and orifice aperture 406 are aligned and adapted to hold the jetting tube 432 and transducer 434 such

that the orifice 433 of the jetting tube 432 extends into the orifice aperture 406.

An electrical transducer activation pulse is supplied to the piezo-electric transducer 434 from the jetting head control unit 500 by means of two contact pins 422. A quantity of fluid will be dispensed from the jetting tube for each applied activation pulse. The activation pulse can be produced by a variety of conventional circuits or commercially available units. Therefore a detailed description of such a circuit will not be provided. However, a circuit for producing a series of activation pulses is provided in the description of the printing embodiment below. Due to the differing constraints involved in dispensing and printing, the circuit in the printing embodiment is not required to produce only a single pulse. However, one skilled in the art could, if desired, modify the circuit to produce a single pulse on demand for use in the dispensing embodiment.

Each contact pin 422 defines an enlarged head 423 which is adapted to contact the respective first and second electrodes 437, 436 located on the outer surface of the transducer 434. Two contact pin holders 414, integral with the housing 402, 404, are positioned to hold the respective contact pins 422 under the pin heads 423 such that each pin head 423 electrically engages the appropriate electrode 437, 436 of the transducer 434. Two contact pin engaging posts 420 extend from the housing 402, 404 opposite the contact pin holders 414 to engage and hold the contact pins 422 against the contact pin holders 414. The ends of the contact pins 422 opposite the pin heads 423 extend through the housing 402, 404 by means of contact pin apertures 421. Since the housing sections 402, 404 are formed symmetrically to one another, the contact pins 422 may be optionally attached above the transducer 434.

In operation, the reservoir 200 containing reagent fluid is fastened to the jetting head 400 such that the fluid supply tube 430 extends into the reagent fluid. The filter 300 may be fitted to the free end of the supply tube 430 or positioned inside the reservoir 200. Air is supplied through the channel 410 around the supply tube 430 to prevent the reservoir 200 from falling below atmospheric pressure. The air is prevented from entering around the supply tube 430 and into the transducer chamber 403 by the seal created between the sealing teeth 412 and the supply tube 430. The jetting tube 432 may be primed by slightly pressurizing the reservoir 200 to cause the reagent fluid to travel through the supply tube 430 and into the jetting tube 432. Once primed, the fluid is prevented from substantially withdrawing from the jetting tube 432 by the surface tension of the reagent fluid at the orifice 433.

The transducer activation pulse is conducted to the contact pins 422 of the jetting head 400. The contact pins 422 communicate the high voltage pulse to the electrodes 437, 436 of the transducer 434 with polarity such that the concentrically mounted transducer 434 expands. The rate of expansion is controlled by the rise time of the high voltage pulse which is preset to generate a rapid expansion. The expansion of the transducer 434 causes the jetting tube 432, which is epoxied to the transducer 434, to also expand. The expansion of the tube 432 generates an acoustic expansion wave interior to the tube 432 which travels axially towards the orifice 433 and towards the fluid receiving aperture 431. When the expansion wave reaches the orifice 433, the reagent fluid is partially drawn inwardly. However, the surface tension of the fluid acts to inhibit substantial inward fluid movement.

When the expansion wave reaches the end 431 of the tube 432, the expansion wave is reflected and becomes a compression wave which travels towards the center of the piezo-electric tube 434. The high voltage pulse width is adapted such that when the reflected compression wave is beneath the piezo-electric tube 434, the high voltage pulse falls, resulting in a de-expansion of the transducer 434 and the jetting tube 432. This action adds to the existing acoustic compression wave in the interior of the jetting tube 432. The enhanced compression wave travels toward the ori fice causing reagent fluid to be dispensed from the tube 432. The fluid is propelled from the orifice 433 as a small droplet 2 and deposited in the selected mixing cell 904 positioned by the transportation unit 902. One droplet 2 is dispensed for each transducer activation pulse. This mode of dispensing is referred to as the drop on demand mode.

In some instances, the droplet 2 may be accompanied by at least one smaller satelite droplet. However, even if satelite droplets are present, the volume and velocity of the reagent droplets 2 are highly reproduceable. This reproduceability allows for precise dispensing of uniform, controllably sized droplets 2 of reagent fluid into the mixing cell 904.

The droplets 2 of reagents impact the mixing cell 904 with sufficient force and volume to cause fluidic mixing of the reagents. Once the desired amounts of the selected reagents are deposited in the selected mixing cell 904, mixing cell 904 is transported to the detection station 906 where the mixed reagents may be extracted for use or analyzed for assay results.

The dispensing system 30 provides numerous advantages based upon the ability of the reagent jetting head 400 to rapidly and reproduceably eject uniform quantities of a wide range of reagents. The reaction times of some chemical processes are dependent upon the volume of the reagents used. The ability of the dispensing system 30 to dispense such minute amounts of reagents thereby reduces the processing time

of certain chemical assays. Furthermore, some chemical assays require a wide rampe of dilution ratios. Many conventional dispensing systems are unable to dispense the reagents in volume small enough to make the desired assay practical. The dispensing system of the present invention overcomes this disadvantage.

In addition to dispensing reagent fluids, certain embodiments may be used for precision printing of reagents onto a printing medium such as filter paper to produce an assay test strip. A printing system 10 using the present invention is represented in Fig. 3. Structure similar in form and function to structure described above will be designated by like reference numerals. The printing system 10 comprises a reagent fluid reservoir 200, a filter 300, a reagent jetting head 400, a jetting head control unit 500, an interface 600, a computer 700, and an x-y plotter 800.

The x-y plotter 800 is a commercially available pen plotter, mechanically modified in a conventional manner such that the pen is replaced with the jetting head 400. The general operation and structure of the plotter 800 will not be described in detail. The plotter 800 accepts commands from the computer 700 thru a standard RS-232 serial interface contained within the interface unit 600. The plotter 800 processes the commands and produces control signals to drive an x-axis motor (not shown) and a y-axis motor (not shown). The x-axis motor is used to position the jetting head 400 and the y-axis motor is used to position a drum (not shown) to which the printing target 1 is attached.

The plotter 800 produces a pen down signal PENDN. This signal is applied to the control unit 500 and indicates that the plotter 800 is ready to begin a printing operation.

The control unit 500 also receives control signals from the interface unit 600. These signals include signals HIGHER', LOWER' to control the magnitude of the pulse applied to the transducer 434; a reset signal RST to reset the control unit 500; and a series of print signals PRT. The generation of these signals will not be described in detail since their production is performed by the conventional interface unit 600.

The jetting head 400 and fluid supply system 200, 300 are initialized and operate substantially as described above. The jetting head control unit 500, shown in Figs. 5a - 5e comprises a print control circuit 510, a pulse generator 530, a high voltage supply 540, and a strobe pulse generator 560. The control unit 500 also comprises a power supply. However, since the power supply is of conventional design it will not be shown or described in detail.

The print control circuit 510 receives the pen down signal PENDN from the plotter 800 and comprises a transistor Q100, a one-shot circuit U100, two NAND-gates U101, U102, a line decoder multiplexer U107 and four Inverters U103-U106. The pen down signal PENDN is applied to the base of the transistor Q100 by resistors R100, R101 and diode D100. The emitter of transistor Q100 is tied to ground and the collector is connected to the +5 volt supply by resistor R102.

The one-shot U100 comprises inputs A, B and an output Q. The B input of the one-shot U100 is connected to the collector of the transistor Q100 and the A input is tied to ground. The time period of the pulse produced by the one-shot U100 is determined by a resistor R104, a variable resistor R105 and a capacitor C100. The output Q of the one-shot U100 is combined with the collector output of the transistor Q100 by the NAND-gate U101 and then inverted by the NAND-gate U102. The circuit is operative to produce an adjustable delay in the application of the pen down signal PENDN to the control unit 500.

The line decoder U107 is circuited to function as a 3 input AND-gate. The output of the NAND-gate U102 is applied to the first input of the decoder U107; the print signal line PRT comprising a series of pulses from the interface unit 600 is applied to the second input; and a jetting head ON/OFF signal from switch S1 is applied to the third input. The inverter U106 inverts the output of the line decoder U107 to generate the print control signal PRT and the inverters U103-U105 invert the control signals LOWER. HIGHER, and RST signals, respectively.

The high voltage supply 540, shown in Fig. 5b, provides + 175 volts DC to produce a maximum pulse of + 150 volts peak to peak at the reagent jetting head 400. The high voltage supply 540 comprises differential amplifier U12 and transistors Q1, Q2, Q13, Q14. A stable reference voltage of -2.5 volts DC is produced at the junction of a reservoir R13, connected to the -15 volt supply, and a diode CR6, connected to ground. The reference voltage is combined with a resistor R14 to produce an adjustable, stable voltage reference for the amplifier U12. The reference voltage is applied to the inverting input of the amplifier U12 through a resistor R11. The noninverting input of the amplifier U12 is connected to ground by a resistor R12. The amplifier U12, in combination with a feedback resistor R10, produces an output signal proportional to the difference of the voltage reference signal and the ground potential.

The output of the amplifier U12 is applied to the base of the transistor Q2 whose collector is connected to the +15 volt supply. The signal produced at the emitter of the transistor Q2 is applied to the base of the transistor Q1 through resistors R8. R6. R5. a transformer L1 and diodes CR4. CR2. CR1. The emitter of the transistor Q1 is connected to ground and the collector is connected to the +15 voltage supply through the

transformer L1. A diode CR3 connects the collector of the transistor Q1 to the junction of the resistor R5 and the diode CR4. The transistor Q1 is biased for proper operation by resistors R7, R6, R5. The resistor R7 and a capacitor C22 connect the junction of the resistor R8, R6 to the +15 voltage supply.

The transistor Q1 and the transformer L1 form a "flyback" blocking oscillator. Any increase in current supplied by the transistor Q1 produces an increase in energy transferred through the secondary winding of the transformer L1 and diode CR5. Therefore, an increase in current supplied by the transistor Q1 results in an increase in power available to the high voltage output. The diodes CR1-CR4 form a "Baker clamp" which prevents transistor Q1 from saturating. The clamp thereby avoids transistor storage time.

The diode CR5 is connected to a multiple pl filter formed by the inductors L3, L2, capacitors C24, C21, C41 and resistors R29. The multiple pl filter attenuates ripple and switching spikes in the signal supplied to the transistor Q13 which produces the high voltage output V++. A resistor R64 connects the base of the transistor Q13 to the emitter and to the resistor U29. The base is also connected to the collector of the transistor Q14 by a resistor R65. The base of the transistor Q14 is connected to the +15 volt supply by a resistor R67 and to ground by a resistor R66. The emitter of the transistor Q13 provides a signal HV SENSE which is fed back to the inverting input of the amplifier U12 through a resistor R9. The high voltage output V++ is produced at the collector of the transistor Q13. The proper biasing of the transistor Q13 is provided by resistor R64 and the biasing circuit comprising the transistor Q14, resistors R67, R66, R65.

The pulse generator 530, shown in Figs. 5d, 5e, comprises an opto-isolator U18, a one-shot U23, a digital to analog (D/A) converter U30 and two binary counters U24, U25. The pulse generator 530 accepts control signals PRT, LOWER', HIGHER', RST and produces the activation pulse which is applied to the transducer 434. In normal operation, the PRT' control signal is supplied to the opto-isolator U18 by a jumper JMP between contact points E5, E6. The opto-isolator U18 is of conventional design and comprises a light emitting diode (LED) circuit and a photo-element circuit. A resistor R15 operates as the load resistor for the LED circuit of the isolator and a capacitor C25 suppresses transient noise on the voltage supply to the isolator U18. The output of the isolator U18 is applied to one input of the one-shot U23 whose time constant is adjustably determined by resistors R38. R25 and a capacitor C30. The pulse from the non-inverting output of the one-shot U23 is fed to the base of a transistor Q9. A resistor R39 sets the approximate base current of the transistor Q9 which is used as a level shifter for converting the CMOS signal level to the +15 volt DC signal level.

The control of the rise and fall rates of the pulse generator 530 is accomplished by directing a pair of current source transistors Q11, Q12 to charge and discharge a capacitor C57. The transistor Q11 is operative as a source of current and the transistor Q12 is operative as a sink for current. A transistor Q10 controls the level of the current by applying an appropriate bias current through a resistor R56 to the base of the transistor Q11. The biasing of the transistors Q11, Q12 is critical to the proper rise and fall rates. Therefore precision voltage references CR13, CR15 are used to provide respective bias reference voltages. A temperature compensation network is formed from zener diodes CR14, CR16 and resistors R55, R54 to maintain stable operation of the transistors Q11, Q12, respectively. The variable resistors R49, R52 may be used to adjust the fall time and rise time, respectively, of the output pulse applied to the reagent jetting head 400. A plurality of resistors R45, R46, R47, R48, R49, R51, R52, R53, R56, R57, R58 are used to properly bias the transistor Q10, Q11, Q12 and capacitors C55, C60 are circuited to maintain stability of the circuit.

The impedance of the output stage of the rise and fall circuitry Q10, Q11, Q12 is very high. With such a high impedance, circuit elements attached to the capacitor C57 could affect the linearity of the rise and fall time constants. Therefore, an FET input operational amplifier U32 is used as an impedance interface. The amplifier U32 is configured in the noninverting mode and circuited with capacitors C58, C59 for stability.

The output of the amplifier U32 is applied to an inverting amplifier U31 by means of a resistor R62. The amplifier U31 inverts and conditions the pulse control signal with the aid of resistors R59, R60. Resistors R61, R63, connected to the -15 voltage supply, provide a means for adjusting the DC level offset of the amplifier U31 output signal. Capacitors C51, C52 are connected to enhance the performance and stability of the circuit.

The output of the amplifier U31 is applied by means of a resistor R41 to the positive voltage reference signal input REF(+) of the D/A converter U30. The negative voltage reference signal input REF(-) is tied to ground by a resistor R40. The D/A converter U30 produces output signals IOUT, IOUT which are proportional to the difference between the positive and negative voltage reference signal inputs REF(+). REF(-). Capacitors C48, C49, C50 are connected to the D'A converter U30 to enhance stability.

The D/A converter outputs IOUT, IOUT are also proportional to an 8-bit binary value applied to inputs B1-B8. The binary value is supplied by the counters U24. U25 which are controlled by the function signals LOWER'. HIGHER' and RST. The LOWER' signal and the HIGHER' signals are applied to the count up and

count down inputs CU, CD of the counter U24 by means of opto-isolators U18, U20, The carry and borrow outputs CY, BR of the counter U24 are connected with the count up and count down inputs CU, CD of the counter U25. The reset inputs RST of both counters U24, U25 receive the RST signal by means of an opto-isolator U21, Resistors R16, R17, R18 are used as load resistors for the LED circuits of the isolators U19, U20, U21 and capacitors C26, C27, C28 are used to enhance the stability of the isolator circuits.

The counters U24, U25 may optionally be preloaded to the selected 8-bit binary value through input lines TP0-TP7. The input lines TP0-TP7 are normally biased to the logical high signal state by resistive network U22. The selected binary value is loaded into the counters U24, U25 by pulling the respective inputs TP0-TP7 low and applying an external, active low, load signal EXT LOAD to pin TP8. The load signal pin TP8 is connected to the load inputs LOAD of the counters U24, U25 and conditioned by a clipping circuit comprised of diodes CR9, CR10 and a pull-up resistor of the resistor network U22.

The noninverted and the inverted outputs IOUT, IOUT are connected to the inverting and noninverting inputs of a differential amplifier U29. The output of the amplifier U29 is fed back to the inverting input by a resistor R50. The amplifier U29 converts the current output of the D/A converter U30 to a voltage output. Capacitors C56, C47 are provided to enhance circuit stability.

The output of the amplifier U29 is applied to the noninverting input of the amplifier U28. The output of the amplifier U28 is fed back to the inverting input by means of a capacitor C46 and a resistor R37. The inverting input is also connected to ground by a resistor R36. To enhance the frequency response of the amplifier U28, a resistor R43 and a capacitor C54 are connected between the frequency compensation input FC and ground. An adjustable DC offset is provided by connecting the output offset inputs OF, OF with a variable resistor R42. The wiper of the resistor R42 is connected to the high voltage power supply output V++.

The output of the amplifier U28 is also connected to the base of a transistor Q4 and through diodes CR11, CR12 to the base of a transistor Q7. The transistor Q4, Q7, Q3 and resistors R30-R35 form an output circuit capable of driving high capacitive loads at high slew rates and wide bandwidth. The variable resistor R31 may be used to set the maximum current through the bias network R30, R33 by measuring the voltage drop across resistor R35.

The strobe generator 560 produces a strobe pulse and comprises transistors Q101-Q105 and a oneshot circuit U108. The strobe intensity is determined by the circuit comprising the transistors Q101-Q104 and resistors R109-R115. The circuit is connected to the anode of the LED 900 and receives two inputs from the interface unit 600 to produce four levels of light-intensity in the LED 900.

The activation aand duration of activation of the LED 900 is determined by the one-shot U108 and the transistor Q105. The one-shot U108 comprises inputs A, B and an output Q. The strobe signal STROBE is applied to the B input from the interface unit 600. The duration of the one-shot U108 output pulse is controlled by the adjustable RC network R107, R108. The output Q is applied to the base of the transistor Q105 by resistor R108. The collector of the transistor Q105 is connected to the cathode of the LED 900 to draw current through the LED 900.

The computer 700, control unit 500 and plotter 800 must be initialized. The initialization of the computer 700 and the plotter 800 will not be discussed since these units are of conventional design and operation.

To initialize the jetting head control unit 500, the computer 700 directs the interface unit 600 to issue a reset command. The reset signal RST is conducted to the control unit 500 whereupon the counters U24, U25 are cleared. The computer 700 then retrieves from its memory, or by conventional operator input, the desired digital setting for the D/A converter. This setting may also be calculated from data and may be tailored to specific sizes of jetting heads 400 or reagent fluids. The computer 700 then issues a series of commands, through the interface unit 600, to increment or decrement the counters U24, U25 to correspond to the desired binary setting. If the command directs that the counters are to be raised, then the HIGHER' signal is applied through the opto-isolator U20 to the count up CU input of the counter U24. Similarly, if the command directs that the counters are to be lowered then the LOWER' signal is applied through the opto-isolator U19 to the count down CD input of the counter U24. Since the carry and borrow outputs CY, BR of the counter U24 are connected to the count up and count down inputs CU, CD, respectively, of the counter U25, the digital setting applied to the D/A converter U30 may range from 0 to 255. Alternately, the counters U24, U25 could be initialized to a desired setting by loading the binary value on the lines TP0-TP7 and strobing the EXT LOAD line.

Once the control unit 500 and the plotter 800 are initialized, the printing cycle may begin. The computer 700 issues a command to the interface unit 600 to produce the series of PRT signal pulses. The computer 700 then commands the plotter 800 to print, for example, a line along a selected path. The plotter 800 positions the jetting head 400 and target 1 and issues the pen down signal PENDN. The signal is delayed by the print control circuit 510 to ensure that the target 1 is properly positioned. At the expiration of the

delay, the signal is ANDed with the closed enable switch \$1 and the series of print pulses PRT. The result of the AND operation is the application of the PRT pulses to the pulse generator circuit 530.

The PRT signal is applied through the jumper JMP to the opto-isolator U18 and then to the one-shot U23. The one-shot U23 produces a pulse signal which is then converted from CMOS signal levels to the 15 volt DC signal level by the transistor Q9. The rise and fall circuitry comprising Q10, Q11, Q12 converts the square wave pulse into a pulse having the rise and fall characteristics preset by the resistors R49, R52. The conditioned pulse is then amplified by the amplifier U32 and applied to the amplifier U31.

The amplifier U31 converts the polarity of the conditioned pulse to that acceptable by the D/A converter U30 and supplies an adjustable DC offset. The DC offset is used to counteract possible distortion attributable to the amplifier U31. The distortion arises in that, for the amplifier U31 to be adequately responsive, a small degree of current must flow through the resistor R41. This current creates an offset condition at the output of the amplifier U29 which is then scaled by the D/A converter U30 in correspondence with the binary data. The resistor R63 allows a small amount of current to be applied to the amplifier U31 to control the offset voltage attributable to the current flowing through the resistor R41.

The D/A converter U30 scales the difference between the inputs REF(+), REF(-) using the binary data supplied to input lines B1-B8 to produce a current output pulse IOUT and a current inverted output pulse IOUT. The two outputs IOUT, IOUT are fed to the amplifier U29 which convert the current outputs into a single voltage output. The scaled, conditioned pulse is then applied to the output circuit comprising the amplifier U28 and the transistors Q3, Q4, Q5, Q6, Q7. The circuit produces a high voltage pulse with the aforementioned rise and fall characteristics to drive the piezo-electric transducer 434.

The high voltage pulse is applied to the transducer 434 and causes a droplet 2 of fluid to be propelled onto the target 1. Since the pen down signal PENDN is still applied, additional droplets 2 are produced from the jetting head 400. The plotter 800 moves the jetting head 400 and target 1 along the desired path during the emission of the droplets 2 to produce the desired printed line. When the printing is complete, the plotter 800 removes the pen down signal PENDN and the droplet emission stops. Of course it should be understood that dots, circles and the like could be produced by appropriate positioning of the target 1 and jetting head 400.

The size and uniformity of the droplets 2, as well as the presence of any satelite droplets, may be observed with the aid of the scope 950 and the LED 900. The scope 950 and the LED 900 are positioned such that the droplets 2 pass between the scope 950 and the LED 900 and within the focal range of the scope 950. The strobe pulse when applied to the LED 900 causes the LED 900 to momentarily flash. The timing of the activation and the width of the pulse may be adjusted such that the flash occurs when the fluid, expelled in response to the high voltage pulse, is between the scope 950 and the LED 900. The dispensed quantity of fluid may then be observed in flight or at or near the moment of separation from the orifice 433. Corrections based on the observation may then be made to the system 10.

Since each droplet 2 is small in volume, the droplet 2 may be rapidly absorbed by the target 1, thereby allowing rapid and precise placement of a variety of reagents on the target 1 with reduced drying time and reduced potential of fluidity mixing. In addition, the ability to place small droplets 2 in a precise manner enables the target 1 to be printed in a high density matrix with a variety of reagents as isolated matrix elements.

In some printing applications, particularly when printing fluids of flow viscosity and surface tension, it may be desirable to force the fluid through the jetting tube 432 under pressure and allow the vibrations produced by the transducer 434 to break the emitted fluid stream into precise droplets 2. Under this mode of printing, the emission of droplets 2 can not be stopped by cessation of the tranducers activation bulse. It is therefore necessary to prevent fluid emission by other means. One preferred means of momentarily stopping emission of the droplets is shown schem atically in Fig. 4. In this arrangement, structure similar to structure represented in Fig. 3 in form and function, is represented by like reference numerals.

The arrangement, generally represented by the numeral 20, includes a closed reagent recirculation system comprising a normally close three way valve 970, a sump 960 and a recirculation pump 980. In the continuous mode, the reagent fluid is forced out the orifice 433 by hydraulic pressure and broken into a series of substantially uniform droplets 2 by movement of the transducer 434. A regulated, filtered air supply 100 is used to pressurize the reagent fluid reservoir 200. The reagent fluid within the reservoir 200 may optionally be agitated by a magnetic stirer unit 990. This is especially useful for reagent fluids comprising suspended particles.

The three-way valve 970 comprises a common channel, a normally open channel and a normally closed channel. The fluid is forced through the filter 300 and applied to the normally closed channel of the valve 970. When the normally closed channel is closed, the normally open channel of the valve 970 functions as a vent for the reagent jetting head 400. The common channel is connected to the reagent supply tube 430

of the jetting head 400. The reapent supply tube 430 is also connected to the samp 960.

20

In operation, the normally closed channel is opened by an appropriate signal supplied by the computer 700 which also closes the normally open channel. When the normally closed channel is opened, fluid is permitted to pass to the sump 960 and to the jetting head 400. The sump 960 collects the reagent fluid not transferred to the jetting head 400. The sump 960 supplies the collected fluid to the inlet side of the recirculating pump 980 which returns the fluid to the reservoir 200. The returned fluid is then mixed with the contents of the reservoir 200 and is available for recirculation.

When operating in the continuous mode, rather than interrupt the continuous stream of print pulses to the jetting head 400, the printing may be momentarily stopped by closing the normally closed channel of the valve 970. The closing of the normally closed channel stops the flow of reagent fluid to the jetting head 400 and allows the jetting head 400 to vent to atmospheric pressure. With the fluid supply blocked, the transducer 434 is unable to expel further droplets 2. Thus, if positioning of the target 1 by the plotter 800 requires a longer time interval than the time between droplet 2 emission, the computer 700 may close the normally closed channel of the valve 970. The plotter 800 may then position the target 1 or position a new target 1 as desired.

When printing, the active ingredient of the reagent is tailored to achieve a desired concentration per unit area on the target 1. However, to a certain extent the final concentration per unit area can be adjusted by varying the density of the droplets 2 printed on the target 1. The preferred embodiment is particularly well suited to this application due to its ability to print precise, discrete pels of reagent.

A second preferred embodiment of the jetting head is illustrated in Figs. 6a-6b and is generally represented as 400'. The jetting head 400' comprises housing formed into three sections 401', 402', 403'. The housing section 403' comprises a recessed region which forms the reagent fluid reservoir 200' when the housing section 403' is positioned against housing section 402'.

The jetting head 400' further comprises a piezo-electric transducer 434' and a reagent jetting tube 432' similar to those of the first embodiment. The jetting head 400' and the transducer 434' are most clearly shown in Fig. 6b. The jetting tube 432' defines an orifice 433' at one end and a reagent fluid receiving aperture 431' at the other end. The transducer 434' is mounted to the jetting tube 432' concentrically about the mid-region of the tube 432' with epoxy.

The transducer 434' and the jetting tube 432' are positioned in channels 420', 418', 416' located in the housing sections 402', 401'. The channel 416' comprises a plurality of sealing teeth 412' operative to engage and seal against the fluid receiving end 431' of the jetting tube 432'. The channel 416' is connected to the reagent fluid supply channel 430'. The supply channel 430' is connected with the fluid reservoir 200' by means of an aperture 431' through the housing section 402', shown in Fig. 6b.

The reservoir 200' comprises a flexible reservoir lining 201' adapted to contain the reagent fluid. The lining 201' comprises one aperture which is connected to the housing 402' to allow the fluid to pass from the lining 201'. A vent (not shown), located in the housing 403', allows the space between the reservoir 200' and the lining 201' to be vented or pressurized. A filter 300' is positioned within the aperture 202' to trap unwanted particulate foreign matter.

Electrical pulses are supplied to the transducer 434' by means of two contact pins 422'. The pins 422' are inserted through respective apertures 419' of the housing section 402' and respective apertures 421' of the housing section 403'. Two thin electrically conductive strips 410', 411', shown in Fig. 6b, are used to connect the transducer 434' with the contact pins 422'. A protective shield 405' extends from the housing position 403' to partially isolate the protruding portions of the contact pins 422'.

The function and operation of the jetting head 400' is similar to that of the jetting head 400 and therefore will not be discussed in detail. The collapsible inner lining 201' of the reservoir 200 allows the jetting tube 432' to be primed by pressurizing the reservoir 200' through the vent 205'. Once primed, the jetting head 400' may be used as described above in reference to the jetting head 400.

The jetting head 400 provides an advantage in that the entire fluidic system is contained in one housing. Such containment allows for fast and efficient replacement of the jetting heads without fluid contamination problems.

A third preferred embodiment of the jetting head is shown in Fig. 7 and generally represented as 400°. The jetting head 400° comprises a housing 403°, a reagent fluid supply tube 406°, a piezo-electric transducer 434° and an orifice plate 404°. The housing 403° defines a conically shaped fluid chamber 432°. An orifice plate 404°, defining an orifice 433°, is fastened to the housing 403° such that the orifice 433° is focated at or near the apex of the conical fluid chamber 432°.

The fluid feed tube 406° is attached to the housing 403° and defines a supply channel 430°. The supply channel 430° is in fluid communication with the fluid chamber 432° by means of a connecting channel 431°. The base of the fluid chamber 432° is formed by the disc-shaped transducer 434°. The transducer 434° is

held in position by a hold down plate 402 attached to the housing 403. The electrical connections to the transducer 434 are of conventional design and are therefore not shown. The housing 403 further comprises a threaded aperture 406 for mounting the jetting head 400.

The letting head 400° operates in a manner similar to the letting heads described above. However, in this letting head the transducer 434° is normally disk shaped. When the electrical pulse is applied, the transducer 434° bends slightly, thereby altering the volume of the conically shaped jetting chamber 432°. The change in volume of the chamber 432° causes the expulsion of fluid through the orifice 433° and the intake of fluid through the supply channel 430° as described in reference to the jetting head 400.

A fourth preferred embodiment of the jetting head is shown in Fig. 8 and is generally represented as 400°. The jetting head 400° is very similar in form and function to the jetting head 400 and will not be described in detail. The jetting head 400° comprises two symmetrical housing sections. The sections may be connected together by means of apertures 409° and screws, not shown. When assembled, the housing sections 404°, 402° form a T-shaped supply channel 410°.

In operation, the jetting head 400" functions in a manner similar to the jetting head 400. The jetting head 400" is especially suited for use in the continuous mode, but may also be used in the drop on demand mode. In the continuous mode, the fluid is circulated continuously through the supply channel 430" allowing the jetting tube 432" to withdraw as much fluid as required.

By way of illustrating and with no limitations intended the following information is given to further illustrate the above described embodiments. The computer 700 is an IBM Corporation Personal Computer with 640 kbytes of RAM memory. The interface unit 600 is a Burr Brown interface unit model number PC 20001. The plotter 800 is manufactured by Houston Instrument as model number DMP-40. Communication between the plotter 800 and the interface unit 600 is performed through a standard asynchronous serial communication port.

The electrical pulse applied to the jetting head 400 to activate the transducer 434 comprises a rise time of approximately 5 usecs, a fall time of approximately 5 usecs and a pulse width of approximately 35 usecs. When the transducer 434 is operated in the drop on demand mode, the voltage potential of the pulse is 60 volts plus or minus 10 volts and the pulse frequency can be up to 4 khz. When the transducer 434 is operated in the continuous mode, the voltage potential of the pulse is 30 volts plus or minus 10 volts and the pulse frequency can be up to 10 khz.

The jetting tube 432 is manufactured from a pyrex glass tube and measures .027 inches outside diameter and .020 inches inside diameter. The tube is drawn to a closed taper in an electric furnace. The tapered end is then cut and ground to a desired orifice opening of .002 to .004 inches in diameter. The tube is cut to a final length of .945 inches in the case of the dispenser embodiment and ultrasonically cleaned in acetone. After being cleaned and dried the large end of the tube is fire polished. If desired, the orifice end of the tube may receive a coating, such as a hydrophobic polymer, to enhance droplet separation from the tube.


The supply tube 430 is formed from .023 inch inside diameter and .38 inch outside diameter polyethylene tubing produced by Intramedic Corp. as model number #14 170 11B. During assembly, one end of the tubing is stretched over a warm tapered mandrel. The stretched end of the supply tube 430 is then inserted over the large fire polished end of the jetting tube 432. The assembly is then cleaned and baked in a circulating air oven at 50°C. for 10 minutes.

The transducer 434 was purchased from Vernitron of Cleveland. Ohio as model number PZT-5H. The electrodes 437, 436 are compnsed of nickel and are separated from each other on the outer surface of the transducer by approximately .030 inches. The jetting tube 432 is inserted into the cylindrical piezo-electric tube 434 and secured with epoxy manufactured by Epoxy Technology of Bellenca. Massachusetts as model number 301. The epoxy is applied at the junction of the tube 432 and transducer 434 with a synnge. The epoxy flows along the tube 432 inside the transducer 434 by capillary action. The assembly is then baked in a circulating air oven at 65°C, for one hour to cure the epoxy.

The contact pins 422 are secured to one of the housing sections 402, 404 with a drop of epoxy. The transducer jetting tube 434, 432 is placed in the housing such that the orifice end 433 of the tube 432 protrudes approximately .030 inches from the housing 403, 404. A drop of silver epoxy is placed between each contact pin 422 and the transducer 434 to ensure a secure electrical connection. Epoxy is also applied to the junction of the housing 402, 404 and supply tube 430. The other section of the housing 402, 404 is then screwed into place.

The periphery of the housing 402, 404 is sealed with a capillary sealer such as cyclohexanone. Epoxy is then added around each contact pin 422 and around the orifice end 433. The assembly is then baked in a circulating air oven at 65°C, for one hour.

The filter 300 is formed from a polyester mesh with 30 um pores and positioned in a polypropylene

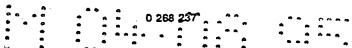

housing. The air pressure supplied to the reservoir 200 during continuous penting operations is regulated at approximately 10 to 30 psi.

The reagents used have the following characteristics:

Printing (drop on demand mode):

Fluid viscosity range: 1 - 30 centipoises
Fluid surface tension: 20 - 70 dyne/cm

Printing (continuous mode):

30

45

55

Fluid viscosity range: up to 50 centipoises

Fluid surface tension: not measured Dispensing (drop on demand mode):

Fluid viscosity range: 2 - 30 centipoises
Fluid surface tension: 20 - 70 dyne/cm

A measure of the performance and selected operating characteristics for a typical jetting head are presented in Figs. 9-11. Fig. 9 is a graph of the mass of a droplet as a function of droplet emission frequency for three fluids. The viscosity of the fluids were 1, 5 and 24 centipoise and the transducer excitation pulse width was 35 microseconds. As shown in Fig. 9, the higher fluid viscosity results in a more stable operating performance of the jetting head. Fig. 10 is a graph of droplet velocity as a function of droplet emission frequency for fluid viscosities of 1, 5 and 24 centipoise. The log of the total fluid weight as a function of the log of the number of droplets emitted is shown in Fig. 11. The fluid used has a viscosity of 2 centipoise, a surface tension of 20 dynes/cm, and a density of .8 grams/cc. The transducer excitation pulse was 80 volts and the excitation frequency was approximately 711 Hz.

Some blood typing reagents and some allergen reagents have very low viscosities and surface tensions. Although in some cases viscosity modifiers, such as glycerol, dextran, glucose, and the like, may be added to increase the viscosity, a few reagents are adversely affected by such modifiers.

Developing stable and reproduceable demand mode jetting is difficult with very low viscosities. Although droplet emission can be established at some fundamental frequencies, the droplets dispensed may have small satelite droplets which reduce the accuracy for metering and dispensing applications. However, even with the satelite drops, sufficient reagent is adequately delivered for most print applications without a substantial decrease in print quality.

Glycerin may be used as a viscosity modifier to improve jetting reliability and to prevent obstruction of the orifice arising from evaporation of the reagent fluid components. Glycerin has been found especially beneficial for those reagents containing particulate material. The evaporation of the fluid component results in a concentration of glycerin located at the orifice. The plug of glycerin substantially prevents further evaporation of the reagent fluid. During the next activation cycle of the transducer, the plug of glycerin is expelled from the orifice.

When operating in the dispensing mode the volume of the droplets can be varied to substantially uniformly contain from 100 pico-liters to 1 micro-liter. The droplets can be produced at a rate of approximately 1 khz to 8 khz. When operating in the printing mode the size of the pel made by each droplet measures approximately .001-.012 inches in diameter.

A copy of the program used in the computer 700 for a printing operation is attached hereto as Appendix A. The values, manufacturer and manufacturing part number of the circuit components of the jetting control unit 500 are substantially as follows:

	Ref. Numeral of Component	Description and Value	Manufacturer and Part No.
	R39,45-48,57,		
	58	RES.10KOHMYWATTS%C.F.	
	R66	RES.1500FMWATTS;;C.F.	
	R3	RES.15KOENWATTS,C.F.	
15	R34	RES.16KOHNWATTS C.F.	•
	R50	RES. 2.4KOHMWATTIM. F.	DALE RL079242G
	R13,23,36,40,		•
	41	RES.2.4KOHMWATTS", C.F.	
	R56	RES. 20KOHN! WATTS .C. F.	
20		• •	
	R8	RES.2200HMWATTSLO.F.	
	R6	RES.270H41WATT5%C.C.	
	R7,12,25	RES_2KOH%WATTSLC.F.	
	R67	RES3.6KOHHWATTS;;C.F.	•
25	R51,53	RES.3.9KOHMWATTS;C.F.	
	R29	RES.300KOHM>WATTS::C.F.	
	R61	RES.30KOYMWATT1;#.F.	DALE RL079303G
	R15-18,26-28,		
	54,55,64	RES.4.7KOHINWATTSUC.F.	
30	R62	RES.45.3KOHMWATTLIM.F.	DALE RN55D4532F
	R30,33	RES. 470H%WATT5%C.F.	
	R21	res.47004Mywatt5;;c.f.	
	R19	RES.47KOHI%WATT5%C.F.	
	R35	RES.5100HM%WATT5%C.F.	
35	R ÷ 3	RES. 6.2KOHMWATT5%C.F.	
	R60	RES.7.5KOHNWATT5%C.F.	
	R37	RES.75KOHMWATT5%C.F.	<u>.</u>
	R9	RES. 76KOHNWATTIXM. F.	DALE RN60D7682F
	R11	RES.8200HNWATTS%C.F.	
40		RES.DIP NETWRK.47KOHM	CT9 761-1R47K
	C21,41,45	CAP.AXIALIME@250VDC	MALLORY #TC56
	C24	CAP.AXIAL220MF@250VDC	MALLORY
		0.	LF2219250C7P3
	ClO	CAP.AKIAL ALUM ELEC.	MALLORY
45	à co	4700 OMF@25VDC	TCG472UO25NIC
	C1,2,3,55,60	CAP. RADIAL DIPPED TANT.	KEMET T350E106M025AS
	053	10MF@25VDC CAP.RADIAL DIPFED TANT.	KEMET
	C53	1MF@35VDC	T350A105K035AS
	C 3.5	CAP.RADIAL DIPPED TANT.	KEMET
50	C36	47MF@10VDC	T350H566MC10AS

Ref. Numeral of Component	Description and Value	Manufacturer and Part No.
C54	CAP.RADIAL SILV MICA	KAHGAN
C57	100PF300VDC CAP.RADIAL SILV MICA	SD5101J301 KAHGAN
¹⁰ C49	20PF300VDC CAP. RADIAL SILV. MICA	SP12200J301 KAHGAN
C39	39FF300VDC CAP.RADIAL X7R MLC	SP12390J301 KEMET
15 C6	.015MF@50VDC CAP.RADIAL X7R MLC	C315C102K1R5CA KEMET C315C223K5R5CA
C30,35,37	.022MF@50VDC CAP.RADIAL Z5U MLC .015MF@50VDC	KEMET C315C153K5R5CA
C4,7		KEMET C315C103K5R5CA
C4,5,6,9,11-19, 22,23,25-28	CAP.RADIAL 25U MLC .22MF@50VDC	KEMET C322C224M5U5CA
C31-34,37,42,43 47,48,50-52	_	
25 C56,58,59	•	
C46 CR7,8,9,10,	CAP.VARI.2-12PF. DIODE SIL.	JOHANSEN #9626 ITT.FAIRCHLD.1N4148
11,12,17 30 CR1,2,3,4 CR5	DIODE SIL.FAST DIODE SIL.FASTHIVOLT	GENL.INST.EGP10D GENL.INST.UF4007
CR6,13,15 CR14,16	DIODE SIL.REF.2,500VDC DIODE SIL.ZENER3.8V.25WATT	NATL.SEMI-LM3852-2.5
U6,13,15,17 ps Q2,9,12	SWITCH 8 POSITION DIP TRANSTOR.COMMON NPN	CTS 206-8 MOTOROLA 2N2222A
08,10,11 Q4	TRANSTOR.COMMON PNP TRANSTOR.HIVOLTHIFREQ.NPN	MOTOROLA 2N2907A . MOTOROLA MPSU10
Q7 Q1	TRANSTOR.HIVOLTHIFREQ.PNP TRANSTOR.HIVOLTHIINPN	MOTOROLA MPSU60 TI, MOTOROLATIP48
40 Q3,14 Q13	TRANSTOR.HIVOLTMPN2N3439 TRANSTOR.HIVOLTPMP	MOTOROLA 2N3439 MOTOROLA MJE5731
U5,27 U23,26	IC 1-SHOT 74HC221 IC 1-SHOT 74LS221 IC COMPARATOR 74HC688	NATL.SEMI MM74HC22IN NATL.SEMI DM741S22IN NATL.SEMI MM74HC688N
U7-10 45 U30 U24,25	IC CONVERTER DAC0800 IC COUNTER 74HC193	NATL.SEMI DACOBOOLCN NATL.SEMI MM74HC193N
U28 U1	IC HI SLEW HI VOLT OF AMP IC HYBRID DC/DC CONVERTER	BURR-BROWN 3584JM BURR-BROWN MODEL 724
U4 50 U3	IC OC DRIVER SN7406 IC OCTAL LATCH 74HC374	NATL.SEMI DM7406N NATL.MM74HC374N
U12,29,31,32 U18,19,20,21	IC OP AMP LF256 IC OPTO ISOLATOR	NATL.SEMI LF256H HEWLTT-PCKRD HCPL2300
R24,42,63 R38,49,52	POT100KOHM%WATT10% POT10KOHM%WATT10%	BOURNS 3622-1-104 BOURNS 3622W-1-103
55 R2O R14,31	POT25KOHM%WATT10% POT2KOHM%WATT10%	BOURNS 3622W-1-253 BOURNS 3622W-1-202

			* * * * * * * * * * * * * * * * * * * *
	Dof Numerical	Description	"ManuTacturer
	Ref. Numeral	and Value	and Part No.
	of Component	and value	
_	VRI	REGULATOR 5VDC	NATL.LM340T-5.0
5	R10	RES. 1MEGOHNWATT5%C.F.	
	R2.4	RES.1.2KOHBWATT5%C.F.	•
	R32	RES.1.6KOHBWATT5%C.F.	
	R44	RES.1.8KOHFWATT5%C.F.	
	R1	RES. 10MEGOHIGWATTS%C.F.	•
10	R5, R22	RES.100HAWATIS%C.F.	
	R65	RES. 100KOHOWATT5%C.F.	•
	R59	RES.10KOHMWATT1%M.F.	DALE RN55D1002F
	R100	RES.2700HI1	
,	R101,108 .	RES.4700HM	•
15	R102,103	RES.1KOHM	.*
	106,109,110		
	R104	RES.47000HM	
	R105	PCT.100KOHM	•
	R107	POT.10KOHM	
20	R111,113	RES.2200HM	
	R112	RES.22OHM	
	R114,115	RES. 470HH	
	C100	CAP.10MF035 VPC	
	C108	CAP.10000 PF	_
25	D100	DIODE	1N4148
	Q100,105	TRANSTOR	2N2222
	Q101,102	TRANSTOR	2N3906
	0103,104	TRANSTOR	2N3904
	U100,U108	IC I-SHOT	74LS123
30	U103,104	IC INVERTOR	74LS04
	105,106		
	U108	IC LINE DECODER	74LS138

Of course, it should be understood that a wide range of changes and modifications can be made to the preferred embodiments described above. For example, the transducer could be of a type other than piezo-electric such as magneto-strictive, electro-strictive, and electro-mechanical. It is therefore intended that the foregoing detailed description be regarded as illustrative rather than limiting, and that it be understood that it is the following claims, including all equivalents, which are intended to define the scope of this invention.

APPENDIX

1 ASE 1

07-14-84 12:24:57

```
Respect Jet Printer
  Respent Calibration
                                                                                                IBN Personal Computer BASIC Commiler V2.00
  Offset bata
                   REB STITLE: Reagent Jet Frinter' SSUSTITLE: Reagent Calibration' SLINESINE: 132
   0030
           6004
    0030
           0004
                    . MITTLE - "REACHL"
    0030
           0004
    4200
           0001
                    "AUTROR - M. A. Enevold
    6030
           0006
           0004
                    COPYRIEST (C) 1985 ABBOTT LABORATORIES
    0020
                    REVISION - 7.0 07-01-84 MAE Microfab modifications
   0030
           8006
                              - 1.0 02-11-84 NAE Creation of initial code
    8638
           6004
    8220
           OCSL
           6006
                    SYSTEM - This code can only be compiled by the BASCON
    0030
                                 EXPILER, It will not run under the INTERPRETER!!
    6029
           5004
    0030
           0004
    0030
                    DESCRIPTION
           0004
                            The reagent calibrate module presents a menu with 12 items arranged
    0020
           6304
                             in 3 columns of 4 rows. The arrow keys allow sovement around the
           0004
    9229
                             table, the + and - keys incresent or decreased values in the first
    0200
           8004
                             column, and the enter key executes commands in the third column.
            8661
    6638
                             The second column is an array of ASCII strings representing respent asses,
    0020
            6664
                             concentration, density, and viscosity. The values entered in column one
    0020
            1000
                             are drop frequency, sulse width, strobe delay, and nozzle musber.
    0030
            6004
                             The cossands in the third column are start/stop, load, save, and exit.
    0030
            8064
    0020
            0005
    0030
            0004
                    'DATA DICTIONARY
                             KENUL
                                           Pointer to which seem item is active (0-11)
    0020
                             NEWS (17,1)
                                            Array for strings used to display the means
    0030
            0004
                             NEXU(17,4)
                                            Array for numbers in the sens display
    6030
            0004
    9020
            0004
                             DIFFI
                                            Differential to save NORUE at arrow key input
                                            Pointer set during warm scan to direct action
            6004
                             TYPEL
    0030
                                            Storage for string imput from mean display
    0030
            6054
                             KEYTUFI
            6004
                             LS
                                            Destination for single taystrole imputs
     0030
                             FILES
                                            String where filensee is built for reagent data file
35
            1000
    0070
                             BERNAMES
                                            String where reagent mase is stored
            6684
     0030
            8004
                                            Row to display special graphics character in semu
    6636
                             Ħ
                                            Column to display special graphics character in sense Special graphics character is read lato here
            6664
                             a
     0030
    0030
            2034
                             17
                             OLD.ARP.VALUEI lateger value for setting pulse amplitude
BIG.VALI Value set to digital port 0 to inc/dec amplitude
     0030
            6454
    0020
            0004
     0030
            0004
                    SOR REASONT. CALIFFATE STATIC
     0030
            0004
     0047
            0004
                             $13. NEXES (17.1) NEXE (17.4)
     0047
            0004
     0048
            OIFE
                             FOSER INITIALIZE:
                                                       'reed imit. values and set screen
    0048
            OIFE
     3400
            OUTE
     BOSE
            DIFE
                             WILE TIPES O S
     0051
            0200
     0051
            0200
                                TIPEL = 8
     0010
                               M = "
    004A
            0204
                                WHILE AS . ""
     00&A
                                  AS . LIXEYS
     0079
            0204
     0083
            0204
                                  IF ACTIVES = 1 AND DOUNTINE ( TIMER THEN GOSUB PEX.DONN
     COAD
            0204
            020A
     0080
```

```
0 268 237
```

5

15

```
20
                                                                                                                                         PAGE 2
07-14-84
   Resgent Jet Printer
Resgent Calibration
                                                                                                                                          12:25:57
                                                                                                     IBM Personal Computer BASIC Compiler V2.00
                    Source Line
   Offset Data
                                                                                   'execute (cr)
25 0080
                                 IF AS = CHRS113) THEN TYPEL = 1:
            020A
                                                                                   'incresent variable
                                 IF AS = "+" THEN TYPEX = 2:
    OOCA
            020A
                                                                                   'decresent variable
                                 IF As * "-" THEN TYPES = 3:
     COEO
            020A
                                                                                   'up arrow key
                                 IF AS = CHRS(D) + CHRS(72) THEN TYPES = 4:
     00F6
            020A
                                                                                   'down arrow key
                                 IF As + CHRS(0) + CHRS(80) THEN TYPEI = 5:
     0113
            020A
                                                                                   'left arm tey
                                 IF AS = CHRS(0) + CHRS(75) THEN TYPEL = 4:
 30 0140
0145
            OZOA
                                IF AS = DKRS(0) + CHRS(77) THEN TYPEI = 7: 'right arrow kay
IF AS > CHRS(47) AND AS ( CHRS(121) THEN TYPEI = 8: 'ascli 0 - 2
            020A
     018A
            020A
     0102
            020A
                                 DH TYPES BOSUB T1, T2, T3, T4, T5, T6, T7, TB
     0102
            020A
     OIDB
            020A
 35 01DF
                               YEX
             070A
                               TYPEL . 0
             020A
     01E6
             OZDA
                               eit sub
            OZOA
     01E6
                      REN SPASE
     OIEA
            020A
```

40

45

50

```
Reagent Jet Printer
                                                                                                                                PAGE 3
   Respect Calibration
                                                                                                                                87-14-66
                                                                                                                                12:24:57
   Offset Data
                                                                                              IEM Personal Computer BASIC Compiler V2.00
                   Source Line
10 01EA
                   ********** SUFROUTINES FOR THIS MODULE *********
           0202
    OLEA
           0204
                                    '(cr) execute consum
    OIEA
           9201
                            IF MEMIN ( 12 THEN TYPES = 0:RETURN: 'exit to print sens, as action
    OIEF
           0204
                            OK NEDUT - 11 SOSUE TIA, TIE, TIC, TIE
    0205
           020C
                            IF REDUCE C IS THEN TYPES . 0
    0214
           6200
    CZZC
           0200
                            BETURN.
    0238
           020C
                                    'start/stop drop flow
                   TIA
    0230
           070C
                            IF NEWDS(12,0) = "START" THEN GOSUS START.INX
IF NEWDS(12,0) = "STOP" THEN GOSUS STOP.INX
    0233
           020C
    025A
           020C
    027F
           0200
                            HERUS (12,0) = TEPS
                            COLOR 0,7:505UB BISPREDEU
    027A
           0210
    OZAC
           0210
                            RETURN
    0280
           0210
                   START. IKK:
    0280
           0210
                            TEMPS . STOP .
    0295
           0210
                                                     'in andule PCI
                            CALL DOT.OM:
    028F
           9710
    0203
           6210
                            LOCATE 17,71:COLOR 27,0:PRINT "PRINTING";
25
                            ACTIVEZ = 1
    02F1
           6210
                            RETURN
    62FB
           6210
    OZFE
           0210
                    STOP. INK:
    OFF
           8210
                            TERPS = "START"
    1020
           0210
                            CALL BOT. OFF:
                                                     'la ecdule PCI
    0303
           0710
                            LOCATE 17,71:COLOR 15,0:FRINT .
    0317
           0210
                            ACTIVEL . 0
    0223
           0210
                            PETILEN.
    0244
           0710
    0348
           0210
                                     "load reagent profile
    0349
            0710
                    TIBE
                             IF REMUSIGATE . . THEN LOCATE 25,1:PRINT "Reagent Hase is not specified";:505UB ANYKEY:RETURN
    0349
            0210
     0371
            0218
            0210
                             EDSUR SEARCH
     0375
            0210
    0317
                             IF II ( (REAMUNE + 1) THEN GOTO FOUND
    0397
            0210
            0214
                             LOCATE 25,10-LEK(MENUS(6,1))/2:PRINT REMUS(6,1); sot Found';
     0343
            0214
                             SUSUB ANYCEY: 'wait for a keybit
    0404
                             RETURN
     0404
            0214
     040E
            0214
     040E
            0214
                    FOUND:
     0413
            0214
                             FILES = RIGHTS (STRS ()1) (LEX(STRS ()1))-1) + "REA, RJP"
                             OPEN FILES FOR IMPUT AS $1:
                                                             'sat pattern data file for read
     0437
            6218
                             INPUT $1,REMICO,01:
                                                     'read frequency
     0448
            0218
                             THPUT #1, NEW (1,6):
                                                      'read applitude
     6443
            0218
                             INCOT $1, RENU(2,0):
                                                      'read strate delay
            6212
     0488
                             1X201 81 XEXU(3,0):
                                                      'read pulse width
     DULE
            0218
                             IMPUT $1,RENU(4,0):
                                                      'read rise time
     0451
            0718
                             1KPUT $1,KENU15,01:
                                                      'read fall time
     04F4
            9718
     0517
            0718
                             1007UT 41,NEXUS (7,1):
                                                      'read concentration
     0519
            0218
                             1XPUT $1,300US(8,1):
                                                      'read density
     0239
            0212
                             14 PUT 41, NDOUS (9,1):
                                                      'read viscosity
      6541
            0719
                             IXPUT $1, NEWS (10,1):
                                                      'read surface tension
     0535
            0218
     OSAT
            0218
 55
```

1版 4

07-14-EL

12:24:57

```
5 Reagest Jet Printer
  Respent Calibration
                                                                                              IBN Personal Computer BASIC Compiler 92.00
   Difset Data
                   Source Line
                                            'done with data file
                            CL 052 111
    0541
           C271
   0220
           0218
                            OPEN "SEADEF.RIP" FOR GUTPUT AS 81
    0280
           0711
                                                             'save filesass is default file
                            PRINT BI,FILES:
           0218
    0203
                                                     'save the directory same as well
                            PRINT BI, MEMIS (6,1):
    0502
           0211
                            LI GOSE III
    05F4
           0218
                                                     'show all paraseters
                            MONE DISP.PARKS:
           0218
    0278
           0218
                            RETURN
   0601
15
           0718
    0102
                            'save reagent profile

IF MENUS(6,1) = " THEM LOCATE 25,1:PRINT "Respent Mase is not specified";:EDSUB ANYMEY:RETURN
                    TIC
           0218
    0102
           0218
    0104
                            OPEN "READIR.RJP" FOR INPUT AS 41
           0218
    CARE
                            INPUT 11,REARINI
    OLSF
           0718
    0471
            0218
                            CLOSE II
                            IF REAMENT ( 80 THEM STOTO SAVE.REA
            0218
    0478
                             LDCATE 25,1:PRINT "Directory is Full (80 reagents saz.)"
            0218
    0187
                             GOSUB ANTICEY: RETURN
            6218
    CAAL
                    SAVE.REA:
    GLAB
            0218
                             SOSUB SEARCH
    0480
            0218
                             IF II > REAMONT THEN GOTO SAVEREAL
25
    4830
            0218
                             REAMONT = 11
     0&C7
            0218
                             COLOR 15,0
     OLCE
            0218
                             LOCATE 25,1:PRIXT NEMUS(4,1); already exists. Replace it with new values? ";
     66DA
            0218
                             AS - "
     070C
            9218
                             WILE AS = ""
            0718
     0716
                                     AS = INCEYS
     0725
            0218
            0218
                             MEXID
     072F
                             LOCATE 25,1:FRINT SPACES (37);
            0212
     073Z
                             IF AS = "Y" OR AS = "Y" THEN GOTO REPLACE
     074F
            0718
                             RETURN
             0213
     0778
     077C
            0218
35
                     SAVEREAL:
             0218
     077C
                             KILL 'READIR.OLD': delete old MARE 'READIR.RJP' AS 'FEEDIR.OLD':
                                                       'delete old backup directory
             0218
     07E1
                                                                       'save old directory
             0218
     0788
                              OPEN "KENDIR. GLD" FOR INPUT AS II
             8218
      0792
                              CPEN "READIR.RJP" FOR OUTPUT AS $2:
                                                                        'set us new dir
             921R
      07A3
      0735
             0218
 40
                                                       'read mumber of dir entries
                              INPUT $1, REAKURE:
      0785
             0218
                              REARUNT - REARUNZ + 1: 'snerease by I
      07C7
             0218
                              WRITE DZ. REAKUNILI
                                                        'save in new directory
             0218
      07E1
             0218
                              FOR I=1 TO REASONS - I
      07E1
             0218
                                  LIKE INPUT 41,AS:
                                                       'read entry from old dir
             0210
      07F&
 45
                                  PRINT $2,A$1
                                                        'write entry in new directory
      0807
             021E
                              METT I
      0817
             0710
      0832
             0220
             0220
      0832
      0825
             8220
                                                       'write see entry to see directory
                               PRINT 12, NEWS (6,1):
      0839
              0220
                                                done with directory
                               CL05E 12:
       6853
              0220
       0842
              0220
       0842
              0220
                               FILES - RIGHTS (STRS (REAMUNE), LEMISTRS (REAMUNE))-11 + "REALRJP"
       0847
              6270
       CORR
```

0220

```
Rageat Jet frinter
                                                                                                                               PAGE 5
  Respont Calibration
                                                                                                                               07-14-64
                                                                                                                               12:24:57
  Offset Data
                                                                                            IBM Personal Computer BASIC Compiler V2.00
10 0223
                           CZEN FILEN FOR OUTPUT AS bin - 'create new pattern data file
          0220
                                                    'store frequency .
   0810
          0220
                           HATTE BI. MEDEU (0.0):
   0013
          0220
                           WRITE BI, MENDILL, BI:
                                                    store applitude
                           MRITE 81, MEDIU(2, 6):
   3080
          0220
                                                    'store strobe delay
                           WRITE 81, MEMU(3, 6):
   OFF
          0220
                                                    'store pulse width
   0915
          8220
                           WRITE BI, MEDBICA, 61:
                                                    'store rise time
15 013F
          0220
                           WRITE 41, NEW (5,0):
                                                    'store fall time
   6440
          8224
   0142
          0220
                           WRITE #1,NEXUS(7,1):
                                                    "store concentration
   9984
          0220
                           STITE $1,RDOUS(E, I):
                                                    "store density
   LAPO
          0220
                           MALTE 01, NEWS (9,1):
                                                    'store viscosity
   0703
          0220
                           WRITE 01, MERCH (10,11)
                                                    'store seriace tension
20 09EA
          0220
   OPEA
                           DLOSE #1:
                                           'done with data file
   09F1
          0220
   09F1
          0220
                           DPEK 'READEF.RJP' FOR OUTPUT AS II
   0402
          0220
                           PRINT $1,FILES:
                                                            'save filesame in default file
   0413
                           PRINT BI, NEWS (4,1):
                                                    'save the directory mame as well
          0770
25 0A33
          0220
                           CLOSE 11
   37.43
          0220
                           RETURN
   0440
          0220
                  SEARCH
   0440
          0220
   OMS
          0220
                           OPER 'READIR.RJP' FOR INPUT AS 81
   CASA
          0220
                           IMPUT BI REAMUNES
                                                    'read number of patterns in dir
30 DASS
          0720
                           11 = 1:
                                                            'set entry pointer
   SALE
          0220
   OAAF
          0220
                  270051
   ,0A74
          0720
                           LINE IMPUT 11,AS:
                                                    'read ment pattern name from dir
                           IF As = MEMUS(6,1) THEN SOTO SEARCH, BOXES
   CARI
          0220
                                                                           'compare name with dir entry
   کلاه
                           11 - 11 + 1
         . 0720
35
                           IF II ( IREAMINE + 1) THEN SOTO SLOOP: check for done
   OME
          0220
   OACI
          0220
                  SEARCH.DOME:
   CACL
                          CLOSE II
          0220
   DACD
                           RETURN
          0220
   CADI
          0220
   OADI
                  713:
          0270
                                   'return with no change to exit reagest calibrate
40
                          PRINT 13, "UH";
   DADE
          0220
   CAEL
          0220
                           CLUSE #31
                                           close cos chansel
   ONED
                          EUZ
          8220
   OFI
          0229
   OUF1
          8728
                  12:
                                   "process "+" key
                           IF MODUL ) 5 THEM RETURN
   OMFA
          0226
45
   6105
                           MENTINE . TIMEN
          0220
   030F
          0224
                           DELTATINE - MENTINE - QUITINE
   OBIF
          0721
                           OUTLINE . REYTLINE
   0921
          OZZ
                           IF DELTATINE > 0.15 THEN HILTS = 1 ELSE HILTS = HILTS + 1
   0141
          02ZE
                           IF MULTS > 100 THEM MULTS = 100
50 0151
          077E
                           MEDELIKERUI,0) = MEDELIKERUI,0) + MEDELIKERUI,3) + MULTI: 'add increseet
                           IF REMINERAL (0) > WENNIMERAL (1) THEN MEMINERAL (0) = MEMINERAL (1): "Check par value
   021F
          622
   0006
          022
                           COLOR 15,1:60SUB BISPADAU:RETURNI
                                                                                    'show wer value
   OCID
          02.TE
   OCID
          077E
                  11:
                                   'process '-' key
   0CZ2
          07<u>2</u>E
                           IF RENUE > 5 THEN RETURN
55 9031
                           RENTINE . TINER
          OZZE
```

```
Respent Jet Printer
                                                                                                                                 PAGE 6
                                                                                                                                 07-14-64
  Respent Calibration
                                                                                                                                 12:24:57
  Offset Data
                                                                                              IBN Personal Computer BASIC Computer V2.00
                   Source Line
                           SELTATINE - RENTINE - OLDTINE
10 0038
          OZZE
                           OLDTINE . RENTINE
   SCAR
          0.7E
                           IF BELTATINE ) 0.15 THEN MULTI = 1 ELSE MILTI = MILTI + 1
   0033
          022E
                            IF MULTI > 100 THEN MULTI = 100
   0C77
          077F
                           REMUMBRUI.0) = REMUMBRUI.0) - REMUMBRUI.TI = MALTI: 'sub incresent
   0027
          077F
                            IF REMU(REPORT, 0) ( REMU(REPORT, 2) THEN REMU(REPORT, 0) = REMU(REMUT, 2);
                                                                                                      'check sis value
   0000
          027E
                           COLOR 15,1:605UB 015PREMU:RETURN:
                                                                                     "show new value
75 0032
          023E
   0049
          022E
                           'process up arrow key

1F MENUI NOD 6 = 0 THEM RETURN:
   0049
          022E
                   14:
                                                                             'is top row already
   SPE
          OZZE
   0063
          022E
                           DIFFI = -1:605UB MEXMENU: RETURN:
                                                                     'sove soister us one
   0074
          0230
20 0074
           0230
                   15:
                                     process down arrow key
                            IF REMUZ KOD 4 = 5 THEN RETURNS
                                                                              Ta bottom row already
   0079
          0230
                            DIFFE . 1:50 SUB NEVRENU: RETURN:
          0230
                                                                              anve pointer foun one
   OBEF
   ODAO
           0230
                                     'process left arrow tay
   ODAO
           0230
                   16:
                            IF INT (NEXUS / &) = 0 THEN RETURN
                                                                     'im left column already
   ODAS
           0230
25 0003
           0230
                            DIFFI . - 6: GOSUB WEYKENU: RETURN:
                                                                     'sove pointer one left
           0230
    ACCO
                   17:
                                     'process right arrow key
    4660
           0730
                            IF INTEREST / 6) = 2 THEN RETURN
                                                                     'is right column already
    ODDR
           0230
                            DIFFI = 4: BOSUB NEVRENU: RETURN:
                                                                              sove pointer one right
    ODFE
           0730
    OEOF
           0230
                                     'input keys into KEYBUFS until (cr) is entered
 30 0E0F
           0320
                   T8:
                            IF KENUE > 10 THEN RETURN
    0E14
           0230
                            LOCATE 25.30:COLOR 31,0:FRINT "ENTER NEW VALUE";:COLOR 15,0
    0E23
           0230
    0235
           0230
                            KEYBUFS = AS
    OESF
           0234
                            WHILE AS () CHRE(13)
                                    LOCATE 25,47:PRINT SPACES(15);
    0E72
           0234
 35 OEBF
           0224
                                    LOCATE 25,47: PRINT TEYBUFS;
   PA30.
           0234
                                    A5 = ""
                                     WHILE AS # **
           0234
    0EB2
                                            AS = INCETS
    OEC2
           0234
                                             IF ACTIVES . I AND DOWNTING C TIMER THEN GOSUB PEN. DOWN
           0234
    OEEE
    CEFA
           0234
 <sup>40</sup> 0EF9
                                     IF AS . CHR. (8) AND LEMIKEYBUFS) > 0 THEN KEYBUFS . LEFTS (KEYBUFS, LEMIKEYBUFS)-1)
           0234
                                     IF AS > CHRS(31) AND LEN(KEYBUFS) ( 15 THEN KEYBUFS - KEYBUFS + AS
    OF 3B
           0234
                            VEXD
           0234
    0F75
           0234
    0F79
                            IF REMUT ) 5 THEN GOTO STORESTRING
           0234
    OF79
    OFB8
           0234
 45 OFBB
                                                     'temp has value of keys input
                            TERP . VALIKEYBUFS)
           0734
    OF98
           023B
                             'round off tesp according to step size in meso array
    OFTE
           0238
                            TERP - INTITERP / (NEMUCRENUL, 31) + .5) + RERUCKEDUI, 3)
    0F98
           023B
    OFDI
           0228
 50 OFD1
                             'test TEMP for easieum and einieum values in emme array
           0228
                            IF TEXP > REMUMERUE. 1) THEN TEXP . REMUMERUE. 1)
           0238
                             IF TENY ( NEW (ADIUL, 2) THEN TENP . MENU (MENUL, 2)
     1019
            0238
     104F
            0738
                             'insert see value into sens array and update screes
     104F
            0238
                            NEWLINERUL.O) . TEMP
     104F
            0238
                            LOCATE 25,30:PRINT SPACES (40):
 55 1043
```

PAGE 7

07-14-44

12:24:57

```
Respect det Printer
   Reagent Calibration
                                                                                                IBN Personal Computer BASIC Compiler V2.00
                             COLOR 0,7:505US BISPREDO
70 1068
            0738
    1014
           0238
                             RETURN
    1072
           0738
    107E
           0232
                             HERLE (HERLT, 1) . KEYPUFE
    1043
            0238
    LOBF
            0238
                             LOCATE 25,30:PRINT SPACES (40);
                             COLOR 8,7: GOSUN DISPREDAL
    1000
            0232
                             RETURN
            0233
    10EE
     10F2
            0233
                    PEL MORE
     10F2
            0228
                             DONOTTINE . TIMER + 1
     10F7
            0733
                             PRINT 13, "3"}
     1107
            0238
                             RETURN
    1117
            0238
     1118
            0238
     1118
            9238
                    ATTETE
                            LOCATE 25,64:PRINT "Strike any key..";
     1120
            0238
                             M = **
     113A
            0238
                             WILLE ME . **
     1144
            0738
                                     AS . INCEYS
    1153
            0738
     1150
            0238
                             LOCATE 25,1:COLOR 15,0:PRINT SPACEF(79);:COLOR 15,1
            0238
     1150
                             RETURN
            873B
     1176
     117A
            0238
                     MEMBERL: 'write ald itee in yellow, point to and highlight new item
     117A
            0238
                             COLOR 14,0:EDSUS DISPHERU
     119F
            0238
                             NEXUL - NEXUL + DIFFI
     1191
            0238
                             IF RENUT . 11 THEN RENUT . 10
     1190
             023B
                              IF REDUL > 15 THEN RENUL . 15
     1105
             0232
                             COLOR 0,7:605UB DISFREND: RETURN
     1181
             0238
     11F7
             0238
     11F7
             0238
                     INITIALITE:
                              'change to second screen and display messages
SCREEM 0,0,1,1:COLOR 7,0:CLS:LOCATE 10,28:PRINT *Imitializing Mean Bisplay';
             0233
     11FC
             0238
                              LOCATE 12, MIPRINT 'Please Wast ....
     1740
             0233
     125A
             0233
                              'Initialize variables
     1254
             0238
      1254
             0232
                              ACTIVES = 01' not printing
      1254
             0233
     1261
             823
                              "imitialize plotter com channel
             0238
      1241
             023
      1261
                              OPEN "COM1:2400,N,B,2" AS #3
     1241
             0232
                              PRINT $3,"1:UECS,EFVI,R"1
             0233
     -1273
             0738
      1293
                               'imitialize digital port
             0738
      1282
                              SERE . 4
             0238
      1283
                              CALL DIGITAL. OUT (SCRE)
      128A
             023A
                              0 = 1822
      127A
             0234
                              CALL DIGITAL OUT (SCREE):
                                                                'pulse reset line to set amplitude to OV.
      1241
             023A
                              5001 = 4
      1231
             023A
                              CALL DISTAL DUT (SCAL)
      1233
             023A
      1203
             OZZA
                               'set bardware sulse width
      1200
             OZZA
                               CALL SET. BOT. MIDTH(5) 'is sodule PCI
      1201
             023A
```

```
Reagont Jet Printer
                                                                                                                                  PASE 1
  Resemt Celibration
                                                                                                                                  67-14-EL
                                                                                                                                  12:24:57
                                                                                               ISA Personal Computer BASIC Compiler V2.00
  Offset- Bata
                  Lower Line
10 12DE
          023
   120E
          OZE
                            'initialize sees arrays
                           RESTORE ARRESTA
   120E
          GZZ
          0235
                           FOR 12=0 TO 17
   1253
                                    READ MERCH (II,0), MERCH (II,1):
   1223
          873Z
                                    READ REPUTET, 1) , REPUTET, 23 , REPUTET, 33 , REPUTET, 41
   1312
          BYTE
                           NETT II
15 137C
          023E
   1325
          023E
                            'set default reagent values
   138F
          923C
   1385
          92X
                           BER (8,0) = 2000;
                                                             'frequency
   1385
          OZX
                                                             asplitute
                           NEW(1,0) = 0:
   1348
          CONT.
                           MENU(2.8) = 1:
                                                             strobe delay
   1324
          EXX
                                                             'pulse width
                           NEDOU (3,0) = 090:
   1320
          023C
                                                             'rise time
                           NEXU(4,0) = 470:
   IJFC
          023C
                                                              'fall time
                           HEXU(5,0) = 070:
   1418
          023
   1429
          623C
   1436
          QZX
                            RESU(4,0) = 0:
                                                             ....
                                                             'concestration
25 1452
          OZX
                           MENU (7,0) = 01
                            EMU(8,0) = 0:
                                                              'deasity
   146E
          023E
                            NEW(7,0) = 0:
                                                              'viscosity
   148A
          DZX
                            班與(10,0) = 0:
                                                                      surface tension
   1446
          0232
   1402
          023E
                            OLD.AXP. VALUES . 0
                                                             'imitial value of 0 volts
   1402
          02X
30
   1409
          OZJE
                            'change active displayed screen to first screen to draw and display parameters
          67<u>5</u>
   1401
   1409
          02¥
                            SCREEN 0,0,0,1:CLS
           073E
   1409
   1486
           07X
                            COLOR 13:LOCATE 1,32:PRINT "REAGENT CALIBRATE";
    1166
           073F
35
                            COLOR 1
           073E
   1507
                            FOR 1-2 TO 79
    150E
           023E
                                    LOCATE 3,1:PRINT "P";:LOCATE 5,1:PRINT "F";:LOCATE 19,1:PRINT "F";
   1518
           07 Y
           623
    134F
                            FOR J=4 TO 18
    1521
           67.X
                                    LCCATE 1,1:PRINT "J";:LOCATE 1,28:PRINT ":";:LOCATE 1,49:PRINT ":";:LOCATE 1,80:PRINT "3";
    1574
           CISE
40
                            MEIT 1
    1508
           02.3E
                            RESTORE TABLE
    1626
           92JE
    1623
           OZSE
                            FOR 1=1 TO 12
                                    READ RI,CI, MILLOCATE RI,CI:PRINT CHRICOTI);
    1137
           123
    1664
           0244
                            MEIT I
    1485
           6244
45
                             'print three headings and instructions
    1485
           0244
           0244
                            COLOR 10.0
    1685
                            LOCATE 4,7:PRINT "DROP PARAMETERS";
LOCATE 4,39:PRINT "REAGENT PARAMETERS"
    1671
           0744
    1643
           0244
                            LOCATE 4,71:FRINT "CONVOOS"
    1603
           0244
    14DF
            0244
50
                             COLOR TILOCATE 21,20:PRINT "Use ";:COLOR 15:PRINT CHR$(27);CHR$(32);CHR$(26);
    1 SOF
            0214
                             PRINT CHRS (37); CHRS (24); CHRS (37); CHRS (25); CCCCR 7: FRINT * to position highlighted cursor*;
            6244
    1729
                             LOCATE 22,18:PRINT "Use "1:COLOR IS:PRINT "+"1:COLOR 7:PRINT " or "1:COLOR IS:PRINT "-"1
    1748
            0744
                             COLOR 7:FRIXT' to scroll current value up or soun";
    178F
            0744
                             LOCATE 23, 25:PRINT "Use "1:COLOR 15:PRINT "DY"1:COLOR 7:PRINT" to activate selection";
    1702
            0244
```

1

20

Respont Jet Printer Respont Calibration. PAGE 9 07-14-84 12:24:57 IBM Personal Computer BASIC Compiler V2.00

	Gffset	Data	Source Line
25			•
	1814	0244	DISP.PARPS:
	1817	0244	'display 18 meau choices in yellow
	1817	0244	
	1317	0244	COLOR 14,0
	1925	0244	FOR MEMUI = 0 TO 17
30	1823	0244	BOSUR DISPRENU
	1831	0244	NEXT MENUL
•	1841	0244	
	1841	0244	'set for reseast mame and highlight it
	. 1841	0244	REDUZ = 6:COLDR 0,7
	1854	0244	ECISIA BIZANIEM
35	185A	6244	
	185A	0244	SCREEN 0,0,0,0
	186F	0744	RETURN
•	1873	0244	REN SPAGE

40

45

50

```
Respent Jet frinter
                                                                                                                                PASE 10
  Respect Calibration
                                                                                                                                07-14-44
                                                                                                                                12:25:57
70 Offset Bate
                   Source Line
                                                                                              138 Personal Computer BASIC Compiler 92.00
   1873
                   LISTEDO:
           2244
                           LEGATE (MEMUL MID 6)+2+7, (INT(MEMUL/6)+28+2)+15+187(MEMUL/12)
    1272
           0214
                           PEINT REMUSERENCE, 0)
   1804
           0244
                            IF NEXUE > 5 THEN GOTO SHOWSTRING:
                                                                  so value to display
    18F2
           0244
75
                           LOCATE (HENUT MID 6)+2+7, MENU (HENUT, 4)
   1901
           9244
                           PRINT USING MEMUS (MEMUZ, 1); MEMU (MEMUL, 0);
    1422
           0244
                           IF REDUCE > 2 THEN RETURN
    1766
           0244
                            ON REDUCE+1 GOSUB SET. FREB, SET. APP, SET. BELAY
    1975
           0244
    4673
           0244
                           RETURN
20 198A
           0244
                   SICASTRINS:
           0244
                           IF MENUI > 10 THEM RETURN
    198F
           0244
                            LUCATE (NEXUZ KOD 61+2+7,48
    197E
           0244
                           PRINT .
    1784
                           LOCATE (HERVI HOD 4)+2+7,48
   1107
           0244
    LYEJ
           0244
                            PRINT MENUS (MENUI, 1)
25 1402
                            RETURN
           0244
    1406
           0244
    1404
           0244
                   SET. FREE:
                            TENP . MENU(0,0)
    HOAL
           0244
                            CALL SET. DOT. MATE (TEMP):
                                                             'is module PCI
    1424
           0244
                            LEDI = 3-INT((TENP+500)/1000)
    1A34
           0244
30 IAST
           0744
                            IF LEDY ( 0 THEN LEDY . 0
                            SCRI = 4 + (LEDI + 32):
                                                                     'set LED intensity
    1869
           0245
                            CALL DISTTAL OUT (SCREEN
           0244
                                                                     'in module PCI
    1489
                            RETURN
    1499
           0246
    IA9D
           0246
    1490
           0744
                   SET. AND:
                            SCRI - CINTINENUINENUI, 0) + 755 / 150):
 35 1AA2
           0246
                                                                             'convert volts to binary number
                            IF SERI = OLD.AMP.VALUEL THEN RETURN
TEMP1 = SERI - OLD.AMP.VALUEL:
    IACB
           0246
    IADC
           0246
                                                                     'calculate delta
                            DLD.AMP. VALUEL . SCRI:
    TAEB
           0248
                                                                     'update old value to current value
                            BIG. VALI = &
    IAEF
           0242
                            IF TEMPE ( O THEN DIG. VALE . 5
    IAFL
           024A
 40 1308
           024A
                            TERPI . ARS(TERPI)
    1915
           0244
                            FOR II * 1 TO TEXPI
    1322
           024C
                                    SCRI . DIE. VALI . (321LEDI)
    133
           0240
                                    CALL DISITAL. DUT (SCRI):
                                                                             'pulse higher or lower
                                    SCRI # 4 + (32 + LEDI)
    134F
           024C
                                    CALL DIEITAL DUTISTREIS
    184F
           024C
                                                                             'sat port to normal
 45 1377
                            EII II
           024C
    1371
           024C
                            RETURN
    1375
           024C
    1895
           024C
                    SET. FELAT:
    139A
           0246
                            TEMP - KENU(2,0)
                            EALL SET.STRONE.DELAY(TEMP):
           024E
                                                             'is sodule PCI
    1334
 50 1164
                            RETURN
           0210
           0245
    IBCA
                    REN SPASE
    1254
           0245
```

PASE II

07-14-64

12:24:57

IBM Personal Computer BASIC Compiler V2.00

```
jest Jet Printer
     Lagent Calibration
   Offset Data
                          ".... THIS HODULE .... BATA USED BY THIS HODULE ....
     LICA
              OZIC
              024C
    IBCA
75 LECA
                          ARRDATA:
              0240
                                                                             Hz","88,888",10000,1,1,18

","888",150,0,1,19

","888",179,0,1,19

","888",779,0,1,19

","888",779,0,1,19
                                     DATA "Frequency
    110
              0740
                                     BATA 'Amplitude
     1351
              0240
                                     DATA "Strake Delay
              9240
     1803
                                     DATA Pulse Width
     1105
              024C
                                    DATA "Rise lise
DATA "Fall lise
BATA "Base", ",0,0,0,0
DATA "Concentration", ",0,0,0,0
DATA "Dessity", ",0,0,0,0
DATA "Viscosity", ",0,0,0,0
DATA "Surface Tension", ",0,0,0,0
DATA "START", ",0,0,0,0
DATA "START", ",0,0,0,0
DATA "LUAD", ",0,0,0,0
DATA "ELIT", ",0,0,0,0
DATA "EXTERNIA", ",0,0,0,0
DATA "T,",0,0,0,0
DATA "T,",0,0,0,0
                                     DATA Rise line
     1307
              024C
20 1909
              024C
              074C
     . 1308
     1200
               024C
     180F
               024C
               0240
     IPEI
               0240
     IBEJ
25 IBES
               024C
      1BE7
               024C
     1869
               024C
      IJEB
               0240
      IBED
               024C
               024C
      HEF
               024C
     19F1
               024E
      1353
               024E
                           TAPLE
      13F3
                                      BATA 3,1,218
      IFF
               024C
      1354
                                      BATA 3,28,210
               024C
                                       DATA 3,49,210
13FC
35 18FE
                024C
                024C
                                       DATA 3,80,191
                0240
                                       DATA 5,1,198
      1000
                                       DATA 5,28,204
                0240
      1002
                                       DATA 5,69,206
                OZ4C
       1004
                                       DATA 5,80,181
                0240
      1004
                                       BATA 17,1,172
       1008
                0245
                                       BATA 17,28,208
                024E
      1COA
                                       DATA 17,57,200
                0240
       1000
                024C
                                       DATA 19,00,217
       ICOE
                0240
       1010
                            END SITS
       1010
                024C
       1017
                024C
 45
       1017
                 0240
                 024C
      50476 Bytes Available
      43960 Bytes Free
```

55

O Warming Error(s) O Severe Error(s)

```
0 268 237
                                                                                         PASE
                 Reagent Jet Printer
                                                                                         07-05-B6
                Pattern Entry/Modification
                                                                                         10:46:13
                                                       IBM Personal Computer BASIC Compiler V2.00
                                 Source Line
                 Offset Data
                                 REM STITLE: 'Reagent Jet Printer' $5UBTITLE: 'Fattern Entry/Modif
                  0020
                         0006
                                 ication'
                                           - *PATENT* Pattern creation, modification, and filing
                                  TENCOK.
                  0030
                         8000
                  0020
                         2000
10
                  0030
                         4000
                                  *AUTHOR - N. A. Enevold
                  0020
                         9000
                                  *COPYRIGHT (C) 1985 ABBOTT LABORATORIES
                         0006
                  0020
                  0020
                          9009
                                  "REVISION - 1.2 03-10-86 WAE Resove Mouse inputs
                  0020
                          4000
15
                                              1.1 02-20-86 RAE Add 80 pattern limit to save
                          9006
                  0030
                                              1.0 01-13-86 MAE Creation of initial code
                          6000
                  0020
                  0030
                          8000
                                           - This code can only be compiled by the BASCOM
                                  SYSTEK
                          8000
                  0020
                                              COMPILER, it will not run under the INTERPRETER!!
                          4000
                  0030
20
                          6000
                  0030
                                  'DESCRIPTION:
                  0020
                          3000
                                          This module allows the user to LOAD, SAVE, DIRectory, D
                  0020
                          0006
                                  RAY and
                                           enter repeat count and other parameters for a pattern t
                   0030
                          0006
25
                                  o be printed.
                                           The low-resolution graphics mode is selected and a menu
                   0030
                          0006
                                   is displayed
                                           across the bottom of the screen. Using arrow keys
                   0030
                          6005
                                           point to the action to be taken and then invoke that ac
                   0020
                          0006
 30
                                   tion with the
                                           Enter key. In the DRAW mode, another menu is
                   0030
                           4000
                                           displayed which allows the user to select from LINE, RE
                          1000
                   0020
                                   Clangle,
                                           Solid RECTangle, or CIRCLE pattern elements.
                   0030
                           0006
 35
                           9009
                   0030
                                   'DATA DICTIONARY
                   0030
                           6006
                                                              51 Row (Elements) by & Column array f
                                           SCHDATZ(50,5)
                           0005
                   0030
                                      storing pattern elements
                                                              Storage for cursor graphics icon
                                           CURSORI (9)
                    0030
                           0001
 40
                                                              Up to 7 senu mases can be saved here
                    0030
                           6000
                                           HENUS (6)
                                           ELNUNZ
                                                              Count of number of elements in a patt
                    0030
                           0006
                                   61
                                                              Current location of graphics cursor
                                           II - YI
                           0006
                    0030
                                                              Value of one dot space on the screen
                                            GRID
                    0030
                           0006
  45
                                    (default is 0.005")
                                                              Location to print instructions
                                            ROYZ COLI
                    0020
                           9007
                                                              Storage for single key-strokes or inp
                           0006
                                            盐
                    0030
                                    at strings
                                                               Which senu is being displayed (1 or 2
                                            RENUKUK
                           9009
                    0020
50
                                                               Pointer to which senu item is highlig
                                            TEM
                           0005
                    0020
                                    hted (0 - 6)
                                                               Mumber of times pattern is to be repe
                                            REPEATZ
                            0006
                     0030
                                    ated when printed
                                                               I and Y axis distance between the pri
  55
                                            10FF YOFF
                     0030
                            0004

    nting of repeated patterns

                                            ROWSP COLSP
                                                               Row and Column spacing for printing a
                     0020
                            6006
```

ultiple sets of patterns

50

55

PAGE 2 Reagent Jet Printer 15 Pattern Entry/Modification 07-05-86 10:46:13 IBM Personal Computer BASIC Compiler V2.00 Offset Data Source Line PATHUMZ 0030 6006 Number of patterns stored in 20 the pattern directory PATDIR.RJP DROWZ DCOLY Row and Column location to display di 0020 6006 rectory entrys NAMES 6000 0030 Pattern mame to be LOADed or SAVEd to directory 25 17 37 0020 0006 Counters used to LOAD or SAVE the ele sent data from/to pattern data file 0030 0004 FILES Wanz of pattern data file 0020 6006 TEPL Which type of element is being drawn. 1 = Line 2 = Rectangle 30 0030 9009 3 = Solid Rectangle 4 = Circle 0030 9009 FLASI Same as TEMPI above STARTMSS: ENDMSS: Message display for startpoint and en 0030 0006 dpoint of element entry 35 Starting cursor position for III YII 0020 0006 element being drawn 0030 9009 DII DYI Delta I and Y values used to re-position cursor after arrow key The bighest number item in th 0020 0006 MATITEM 40 e current senu display 0020 9006 IS IE Starting and ending I position of the menu highlighting blue box 0020 6006 RADIUSI The calculated radius of a ci rcle to be displayed 45 0020 4000 REN SPASE

```
Reagent Jet Printer
                                                                                          PAGE 3
                                                                                          07-05-86
                 Pattern Entry/Modification
                                                                                          10:46:13
                                                       IBM Personal Computer BASIC Compiler V2.00
70
                                  Source Line
                 Offset Data
                                  SUB PATENTRY STATIC
                  0020
                          8000
                  0047
                          0006
                                          WIDTH 40: SCREEN 1:CLS
                  0047
                          0006
15
                  00SF
                          4000
                                          DIM SCADATZ (30,5), CURSORZ (9), MEHUS (6)
                                          ELKURY = 0:17=0:Y1=0:6RID = 0.005
                  0040
                          029A
                          02A4
                  007F
                  007F
                          02A4
                                          LIKE (0,0)-(6,6),,B
                                          LINE (0,3)-(6,3),,B
                          02A4
                   00A1
                                          LIKE (3,0)-(3,4),,B
20
                          0244
                   OOCS
                          02A4
                                          PRESET (3,3)
                   00E9
                          0284
                                          BET (0,0)-(6,6), CURSORI
                   00F3
                          02A4
                   0116
                   0110
                          02A4
                          02A4
                                          LIME (0,0)-(319,190),,B
25
                   011D
                          02A4
                   0140
                   0140
                          02A4
                                           RESTORE INSTRUC
                                           FOR 1=1 TO 4
                   0147
                          0244
                          02A4
                                                   READ ROWI, COLI, AS
                   0151
                                                   LOCATE ROWI, COLI: PRINT AS;
30
                   0164
                          02AC
                          OZAC
                                           NEIT I
                   0180
                   019B
                          0280
                   017B
                                  FIRST:
                          0280
                   01A0
                          0280
                                           KENUKUK = 1
                                           BOSUB SUBMEXU
35
                   OIAA
                          0284
                   01B0
                          02B4
                                           ON ITEM + 1 GOTO PATOIR, PATLOAD, PATSAVE, PATDRAW, REP
                   0180
                          0284
                                   EAT, PATELT
                   OICD
                          .02BB
                                           GOTO FIRST
                   0100
                          0288
 40
                                   REPEAT:
                   CIDO
                          0258
                                                                   'erase blue box around DIR
                   0105
                          0288
                                           GOSUB ITEMBOIERASE:
                                           LOCATE 25,1:PRINT SPACE: (39); 'erase senu line
                          0288
                   OIDB
                                           LOCATE 25,1: INPUT; "Enter Repeat Count ", REPEATI
                   01F8
                          0288
                                           LOCATE 25,1:PRINT SPACES(39); 'erase senu line
                   0218
                          OZBA
 45
                   0235
                          02BA
                                           LOCATE 25,1: INPUT; "Enter I Axis Offset ", IOFF
                                           LOCATE 25,1:PRINT SPACEF(39); 'erase menu line
                          02BE
                   0255
                   0272
                          OZBE
                                           LOCATE 25,1: IMPUT; "Enter Y Axis Offset ", YOFF
                                           6010 FIRST
                          02C2
                   0292
                   0296
                           02C2
                                   PATEITI
50
                                           WIDTH BO:SCREEN O:CLS
                   029B
                           0202
                   0282
                           0202
                                           EXIT SUB
                                   REM SPAGE
                   02B&
                           0202
```

```
Reagent Jet Printer
                                                                                          PASE
                  Pattern Entry/Modification
                                                                                          07-05-86
10
                                                                                          10:46:13
                  Offset Data
                                   Source Line
                                                        IBN Personal Coaputer BASIC Cospiler V2.00
                                                            "list directory of patterns
                   0286
                           0202
                                   PATDIR:
                                                                   'erase blue box around DIR
                   OZBB
                           02CZ
                                           SUSUB ITEMEDIERASE:
15
                                           LOCATE 25,1:PRINT SPACEs(39); 'erase senu line
                   02C1
                           02C2
                                           OPEN "PATOIR.RJP" FOR INPUT AS $1:
                   OZDE
                           02C2
                                                                                   'open directory
                                    file
                   02EF
                           0202
                                           INPUT $1, PATMUMI:
                                                                    'read number of patterns in dir
                                   ectory
20
                    1020
                           0204
                                           LINE (1,1)-(318,189),0,8F:
                                                                            'erase graphics tablet
                    0326
                           0204
                                           1 = 0:
                                                                            'set counter.
                    0220
                           02C4
                    0220
                           0204
                                   DISLOUP:
                                           I = 1 + 1:
                                                                    'set for next value
                    0335
                           0204
25
                    0344
                           02C4
                                           IF I > PATHUNZ THEN GOTO DIREXIT:
                                                                                   'test for done
                    035B
                                           IF INT ((1-1)/44) <> (1-1)/44 THEN GOTO SHOWNEXT
                           0204
                           02C4
                                           IF INT ((I-1)/44) < 1 THEN GOTO SHOWNEXT
                    0384
                    OJAF
                           0204
                                           LOCATE 25,1:PRINT "More to Display. Continue ? (Y or N)
                    03A9
                           0204
30
                                           GOSUB CORLOOP: 'wast for Y or N response
                    03C3
                           02C4
                    0309
                           0204
                                           IF As = "X" THEN GOTO DIRETIT: 'if N then don't contin
                                   ĦĒ
                    OZBC
                           0204
35
                    03DC
                           02C4
                                         · LIKE (1,1)-(318,189),0,8F;
                                                                            'erase graphics tablet
                           0204
                    0401
                    0401
                           02C4 _
                                   SHOWNEXT:
                                            DROWI = ((I - 1) NOD II) + 2: 'calculate row for disp
                    0406
                           0204
                                   lay
 40
                    0422
                           0206
                                            DCOLI = 4:
                                                                            'set column to 4
                    0429
                           02CB
                                            IF ((I - 1) NOD-44) ) 21 THEN DCOLX = 23: reset coluen
                                    if necessary
                    044C
                           0208
                    044C
                                            LINE INPUT $1. AS:
                           0208
                                                                    'read mext name from directory
                                            LOCATE DROWN, DCOLI: PRINT AS;
 45
                    0459
                           0208
                                                                            PRINT NAME
                                            BOTO DISLOOP
                     0475
                           02CB
                     0479
                           0208
                                    DIRELIT:
                     0479
                           0208
                     047E
                            02C8
                                            CLOSE $1:
                                                             'terminate access to PATDIR.RJP
                                            GOTO FIRST
 50
                     0485
                            0208
                     0489
                            0208
                                    REN SPAGE
                     0489
                            02C8
```

```
0 268 237
                                                                                         Proe 5
                 Reagent Jet Printer
                                                                                         07-05-86
                 Pattern Entry/Modification
                                                                                         10:46:13
                                                       IBM Personal Cosputer BASIC Compiler V2.00
                 Offset Bata
                                 Source Line
                  0489
                         0203
                                 FATLDAD:
                                                                  'erase blue box around DIR
                                         BOSUB ITEMBOIERASE:
                         0203
                  OABE
                                         OPEN "PATDIR.RJP" FOR INPUT AS $1
                  0494
                         02CB
                                         INPUT $1.PATKUMI:
                                                                  'read number of patterns in dir
                  04A5
                         0708
                                         FOSUB SETKAME:
                                                                  'prompt for and input patters n
                  0487
                         G2C8
                                 ut
                                          LINE (1,1)-(318,189),0,BF:
                                                                          'erase graphics tablet
                 . O4BD
                         0208
                  04EZ
                         0203
                                          SUSUB SEARCH
                         0208
                  OSEZ
                  04E8
                          02CB
                                          IF II ( (PATHUMI + 1) THEN SOTO FOUND
                         0208
                  DAFR
                                          LOCATE 10,16-(LEN(NAMES)/2):PRINT NAMES; " not Found";
                  O4FC
                          02CA
                                          LOCATE 12,14:PRINT "Strike Any Key"
                  0531
                          02CF
                                          SDSUB ANYKEY: 'wait for a keyhit
                  054B
                          OZCE
                                          SOTO FIRST
                   0551
                          02CE
20
                   0555
                          02DE
                                  FOUND:
                  0555
                          OZCE
                                          FILES = RIGHTS(STRS(IZ), LEN(STRS(IZI)-1) + "PAT.RJP"
                  055A
                          02CE
                                          OPEN FILES FOR INPUT AS $1:
                                                                           'set pattern data file
                   057E
                          0202
                                  for read
25
                                                                   "read number of elements in pat
                                          INPUT $1,ELNUMY:
                   058F
                          0202
                                  tern
                                                                   'read grid size
                          0202
                                           INPUT DI, ERID:
                   05A1
                                                                   'read repeat count
                          0202
                                           INPUT #1, REPEATZ:
                   0583
                                           INPUT $1,10FF:
                                                                    "read x axis offset for repeat
                   0505
                          0202
30
                                                                   'read y axis offset for repeat
                                           INPUT $1, YOFF:
                   0507
                          0252
                          0207
                   OSE9
                                           FOR II = 0 TO ELHUMI - 1
                   OSES
                          0272
                                               FOR JI = 0 TO 5
                   05F7
                          0204
                   OSFD
                          0204
                                                   INPUT $1,30%CATE(II,JI): read file into screen
35
                                   array
                                               NEIT JI
                   0621
                          0206
                                           NEIT II
                   0631
                          0206
                          0206
                                           CLOSE 11:
                                                            'done with data file
                   0643
                   064A
                          0206
40
                                           OPEN "PATDEF.RJP" FOR DUTPUT AS $1
                          0206
                   0649
                                           PRINT #1,FILES:
                                                                          . 'save filename in defau
                   062C
                           0206
                                   It file
                                           PRINT 61, MAMES:
                                                                            'save the directory mas
                   3440
                           0205
                                   e as well
 45
                           0206
                                           CLOSE II
                    067C
                           0206
                    2830
                                            GOTO REDRAW
                           0274
                    0687
                           0208
                    0687
                                   SEARCH:
                    0687
                           02I:6
50
                                                                             'set entry pointer
                                            17 = 1:
                           0206
                    2890
                           0256
                                    SLOOP:
                    0693
                                            LIKE INPUT 41,A5:
                                                                    'read next pattern name from di
                           0226
                    0698
                                            IF AS = MAKES THEN GOTO SEARCH.END:
                                                                                     'compare name w
                    06A5
                           0204
 55
                                    ith dir entry
                                            17 = 17 + 1
                            0236
                     0688
                                            IF II ( (PATKUMI + 1) THEN SOTO SLOOP: check for done
                     0601
                            0708
```

15

0604

0206

SEARCH.END:

15

20

s Reagent Jet Printer
Pattern Entry/Kodification

FASE 6 07-05-86 10:46:13

Offset Data Source Line

IEM Fersonal Computer BASIC Cozpiler V2.00

'not found so close file and display me

0609 0234 ELDSE 81: 55898
0660 0736 RETURN
0664 0236
0664 0206 REN \$PAGE

35

30

40

45

50

55

```
PASE
                  Reagent Jet Printer
                  Pattern Entry/Modification
                                                                                           07-05-86
                                                                                           10:46:13
                                                         ISM Personal Coaputer BASIC Coapiler V2.00
                                   Source Line
                   Offest Data
                                   PATSAVE:
                   06E4
                           0208
                                            SOSUB ITEMSOIERASE: 'erase blue box around DIR
                   OAFQ
                           0706
                                            IF ELMINI . O THEN GOTO FIRST: 'no elements in pattern
                   OSEF
                           0206
                                           OPEN "PATDIR.RJP" FOR INPUT AS $1
                   OFE
                           0206
                                           INDUT #1, PATRUMI
                   070F
                           0206
                                            IF PATHUMI ( 80 THEN BOTD SAVE.PAT:
                                                                                     'directory fall
                           0206
                   0721
                                    at 80 patterns
                                            CLOSE $1
                    0730
                           0206
                                            LOCATE 25,1:PRINT SPACE$(39);:
                                                                                     'erase bottom 1
                           0206
                   0737
                                   ine
                                            LDCATE 25,1:PRINT "Directory is full (80 patterns max)"
                           0206
                    0754
                                   i
                           0206
                                            BOSUB ANYKEY: BOTO FIRST
                    078E
                           0206
                                   SAVE.PAT:
                    077B
                                            SOSUB GETHAME: 'prompt for and get pattern mame
                    0770
                           0206
                    0783
                           0206
                                            GOSUB SEARCH
                                            IF II > PATHUMI THEN SOTO ADD. NEW. PATTERN
                    0789
                           0206
                    079A
                           0204
                                            LINE (1,1)-(318,189),0,8F:
                                                                            'erase graphics tablet
                                            LDCATE 10,13-(LEN(NAMES)/2):PRINT NAMES; already exist
                    07BF
                           0206
                                    5.*;
                                            LOCATE 12,15:PRINT "Replace it?"
                    07F4
                           0206
                                            PATNUMI = 11
                           0206
                    080E
                                            AS = **
                    OB15
                            0206
                                            WHILE AS = ""
                    081F
                            0205
                                                     AS = INKEYS
                    082E
                            0206
                    0828
                            0206
                                            IF As = "Y" CR AS = "y" THEN BOTO SAVE.PATTERN
                            0206
                    0838
                                            BOTO FIRST
                            0206
                    0864
                    8380
                            0206
                                    ADD. NEW. PATTERN:
                    8480
                            0205
                                            KILL "PATDIR.GLD":
                                                                     'delate old backup directory
                     084D
                            0206
                                             NAME "PATOIR.RJP" AS "PATOIR.OLD":
                                                                                      'save old direc
                            0206
                     0874
                            0206
                                             OPEN "PATDIR.OLD" FOR INPUT AS $1
                     OB7E
                                             DPEN "PATGIR.RJP" FCR OUTPUT AS $2:
                     088F
                            0206
                                                                                      'set up new dir
                                             IMPUT $1.PATHUMI:
                                                                      'read number of dir entries
                            0206
                     OBAI
                                             PATRUMI = FATRUMI + 1: 'increase by 1
                     0883
                            02D6
                                             WRITE #2, PATHUKZ:
                                                                      'save in new directory
                     OBBC
                            OZDA
                                             FOR 1=1 TO PATHUME - 1
                     Q2ED
                            0206
                                                 LINE INPUT $1,48:
                                                                      'read entry from old dir
                     08E6
                            OZDA
                                                 PRINT $2,A5:
                                                                      'write entry in new directory
                     08F3
                            02DA
                                             NETT I
                            02DA
                     0903
                                             PRINT $2, NAMES:
                                                                     · 'write new entry to new directo
                            OZDA
                     091E
                                     TY
                                             CLOSE #1:CLOSE #2:
                                                                      'done with directory
                     092E
                             OZDA
                                     SAVE.PATTERH:
                             02DA
                     0930
                                             FILES = RIGHTS(STRS(PATXUMI), LEN(STRS(PATNUMI))-1) + *P
                             02DA
                     0941
                                     AT.RJP*
                                              OPEN FILES FOR GUTPUT AS $1:
                                                                            icreate new pattern dat
                     0945
                             02DA
                                     a file
55
                                              WRITE FI, ELNUMZ:
                                                                      'store number of elements
                      0977
                             02DA
                                                                      'store grid dimension
                                              WRITE BI.GRID:
                      0988
                             OZDA
                                                                       'store repeat count
                             02DA
                                              MRITE 11, REPEATA:
                      0998
                                              WRITE $1,10FF:
                                                                       'store x axis offset for repeat
                      09A9
                             02DA
```

20

25

30

35

40

45



PAGE B Reagent Jet Printer 20 07-05-86 Pattern Entry/Modification 10:46:13 IBM Personal Computer BASIC Compiler V2.00 Source Line Offset Data 'store y axis offset for repeat WRITE BI,YCFF: 0989 0758 **2**5 FOR II = 0 TO ELEUXI - I 0909 025A FOR JI = 0 TO 5 0907 OZDE WRITE \$1,5ENDATZ([Z,JZ): 'write screen a 0220 OPDD rray to file MEIT JI CZEC 0A00 30 MEIT II OALO OZEC CL 35E \$1: 'done with Wata file ·0A22 DZDC OPEN "PATDEF.RJP" FOR OUTPUT AS \$1 0A29 CZDC PRINT \$1,FILES: 'save filename in defau OASB 020C lt file 35 'save the directory name PRINT \$1, NAMES: OA4B 3200 e as well OA5B OZDE CLOSE 11 6010 FIRST 0A62 02DE 0466 OZDE REM SPASE

15

40

45

50

55

```
0 288 237
                 Reagent Jet Printer
                                                                                         PASE
                 Pattern Entry/Modification
                                                                                         07-05-86
                                                                                         10:46:13
                 Offset Data
                                  Source Line
                                                       ISX Personal Computer BASIC Compiler V2.00
5
                         6230
                                  PATERANE
                  0260
                  ESAO
                         STDC
                                          GOSTE TERSOTEFASE
                  0A71
                                          LIKE (1,1)-(318.1891,0,8F:
                          0200
                                                                           Erase graphics tablet
                  0496
                         0270
                  OAPS
                          0200
                                  WEITEL:
10
                  OA?P
                                          MENUNUM = 2
                          OZEC
                   0445
                          0200
                                          ESSUB SUPPEDIU
                  CAAB
                          OZDC
                  CARB
                                          CN ITEM + 1 EDTO ALINE. RECT, SRECT, ACIRCLE, REDRAY, B
                          OZDT
                                  ACKUP
15
                   BEAO
                          OZEC
                                          BOTH KEXTEL
                   CACB
                          0200
                   OACE
                          OZEC
                                  BACKUP
                   OADO
                          OZDE
                                          GOSUB ITEMBOIERASE
                                          GOTO FIRST
                   ÓАŨЬ
                          025C
20
                   OADA
                          0200
                                  ALIKE:
                   OADA
                          OZDC
                   OADF
                          OZDE
                                          TEMPZ = 1
                                          STARTHSES = "STARTING ENOPOINT"
                   ORES
                          OZDE
                   OAFO
                          022
                                          ENDMSES = "ENDING ENDPOINT "
25
                   ORFA
                          0ZES
                                          BOTO EXTERELEMENT
                   OAFE
                          02E4
                   OAFE
                                  RECT:
                          02E&
                   0803
                                          TEMP2 = 2
                          02E6
                                          SOTO RECTASS
                   4080
                          OZEÁ
30
                   OBOE
                          02E6
                                  SAECT:
                   ORDE
                          0256
                   0813
                          OZEA
                                          TEMP2 = 3
                                  RECTASE:
                          02E6
                   OBIA
                   OBIF
                                          STARTINGS = "STARTING CORNER"
                          02E6
35
                                          ENDASSS = "ENDING CORNER .
                   0829
                          OZES
                                          BOTO ENTERELEMENT
                   0823
                          07E4
                   0837
                          03E6
                                  ACTROLE:
                   0837
                          93EP
                   OB3C
                          02E5
                                           TEXPZ = 4
 40
                                           STARTMESS = "CENTER OF CIRCLE"
                   0843
                          02E6
                                           EXCHSES = "POINT ON CIRCLE "
                   OB4D
                          02ES
                   0557
                          02Es
                                   ENTERELEMENT:
                   0E57
                          02E&
                                           SOSUB ITEMBOIERASE
                   OBSE
                           02E6
 45
                                           FLASZ=0
                   0862
                          02Eb
                                           LOCATE 25,1:PRINT SPACES (39);
                   0844
                           OZEB
                                           LOCATE 25,1:PRINT STARTKSSS;
                   0886
                           OZEB
                                           GOSUB DISPCURSOR
                           02EB
                   OBAO
                                   FINDSTART:
                   OBA6
                           02EB
 50
                                           BOSUB XCUSEACT
                   OBAB
                           OZEB
                                           IF AS = CHR$(27) THEN GOTO ABORT
                   OBB:
                           02EB
                                           IF AS = CHR$(13) THEN GOTO SETSTART
                           0288
                   ORCE
                                           GOSUB CURSORMOVE
                   OBOF
                           02EB
                                           GOTO FINDSTART
                   0852
                           OZEB
 55
                    SEE
                           0228
                                   ABORT:
                                           GOSUB FLACECURSOR
                    OBED
                           02E8
                    OBF3
                           OZEB
                                           SOTO HEITEL
```

QBF7

50

55

```
Reagent Jet Printer
                                                                                          PAGE 10
                  Pattern Entry/Modification
                                                                                          07-05-86
                                                                                          10:46:13
15
                  Offset Data
                                   Source Line
                                                        IEM Personal Cosputer BASIC Cospiler V2.00
                    OBF7
                          OZEB
                                   SETSTART:
                                           LOCATE IS, 1: FRINT ENGASES:
                    OBFC
                          0ZEB
                    0016
                           OZEB
                                           FLAGI = TEMPI: III = VI: YII = YI
                    0028
20
                           OZEC
                                           IF FLAGE = 4 THEN PSET (II+4,YZ+4)
                    ೦೮ಽ
                           OZEC
                                   FINDEND:
                    OC5A
                          02EC
                                           BUSUB MOUSEACT
                    0630
                           02EC
                                           IF AS = CHR$ (27) THEN BOTO CANCELEL
                    0C77
                           OZEC
                                           IF AS = CHR$(13) THEN GOTO SAVEEL
25
                    0C8E
                           OZEC
                                           GOSUB CURSORNOVE
                    0094
                           02EC
                                           BOTO FINDEND
                    0097
                                   CANCELEL:
                           OZEC
                    0090
                                           SOSUB PLACECURSOR
                          OZEC
                    OCA2
                                           ON FLASI GOSUB ERI, ER2, ER3, ER4
                           02EC
30
                    OCB2
                           OZEC
                                           FLASZ = 0
                    OCPA
                          02EC
                                           SOTO NEXTEL
                    OCBE
                           02EC
                                   SAVEEL:
                    ಯ
                           02EC
                                           SOSUB PLACECURSOR
                    0009
                           OZEC
                                           IF FLASE = 4 THEN CIRCLE (111+4, Y11+4), 5QR((11-111)^2+(
                                   YI-Y1I)^2},,,,i
35
                    0032
                           OZEC
                                           GOSUB CORRECT
                    0038
                                           IF AS="N" THEN GOTO REDRAY
                          OZEC
                    0D4B
                                   STOREEL:
                           OZEC
                                           SCHDATI(ELHUHI,0) = FLAGI
                    0030
                           OZEC
40
                    ODSA
                           OZEC
                                           SCHDATI(ELNUMI,1) = 111
                    0085
                           OZEC
                                           SCHDATT(ELNUMI,2) = YII
                                           SCHGATI(ELMUMI,3) = II
                    ODAO
                           OZEC
                    ODBB
                                           SCHOATI(ELHUMI,4) = Y:
                           OZEC
                                           SCHDATI(ELMUNI,5) = 5
                    9DD9
                           02EC
                    ODEF
                           OZEC
                                           ELKUNI = ELNUNI + 1
                    ODFB
                          02EC
                                           FLASZ = 0
                    ODFF
                           OZEC
                                           SOTO KETTEL
                    0203
                          OZEC
                                   REM SPAGE
```

```
0 268 237
                 Reagent Jet Printer
                 Pattern Entry/Mocification
                                                                                          07-05-86
                                                                                          10:46:13
                                                        IBN Personal Computer BASIC Commilier V2.00
                 Offset Data
                                  Source Line
                  0503
                          OZEC
                                  REDRAY:
                                          BOSUB ITEMSOTERASE
                  0E0B
                          02EC
                  OEOE
                          02EC
                                          LINE(1,1)-(318,189),0,BF
                                          IF ELKUME . O THEN BOTO NEITEL
                  0E33
                          02EC
                  0E42
                          02EC
                                          FOR 1=0 TO ELNUME-1
                   0E42
                          02EC
                                                   ON SCHDATI(1,0) GOSUB RD1, RD2, RD3, RD4
                  0E5B
                          02F0
                                          KEIT I
                   0E81
                          02F0
                                          BOTO NEITEL
                   OE9C
                          02F0
                   0EA0
                          02F0
                   OEAO
                          02F0
                                   'sesses Sub-routines called by main module sessess
                   0EA0
                          02F0
                   OEAO
                          02F0
                                  SUBMENU:
                   0EA5
                          02F0
                          02F0
                                           LOCATE 25,1:PRINT SPACE$ (39):
                   0EAS
                                           ON MENUNUM GOSUB MENUL, MENUZ
                   CEC2
                          02F0
                   0ED1
                          02F0
                                           FOR I=0 TO &
                   OED1
                          02F0
                                                   READ MENUS (1)
                   OEDB
                          02F0
                   OEF2
                          02F0
                                                   LOCATE 25. (1+6)+2:PRINT MENUS(I);
                                           NEIT I
                   OF2B
                          02FG
                   0F46
                          02F0
                   0F46
                                           READ MAILTEN
                          02F0
                                           ITEM = 0
                   OF4D
                          02F4
                   0F57
                          02F4
30
                   0F37
                          02F4
                                   HEWITEN:
                                           BOSUB HEWITEMBOY
                   OF5C
                          02F4
                   0F62
                          02F4
                   OF62
                          02F4
                                   NEITITEM:
                                           BOSUB ITEMSEARCH
                          02F4
                   0F67
35
                                           IF AS = CHRS(13) THEN RETURN: TIEN has correct value
                   OFAD
                          02F4
                                           IF LEN (AS) C 2 THEN BEEP: BOTO NEXTITEN
                   OFB4
                           02F4
                   OF9A
                           02F4
                                           IF ASCINICS(AS, 2.1)) = 75 THEN BOTO LEFTAR
                                           IF ASCINIDS(AS.2.1)) = 77 THEN BOTO RIGHTAR
                   OFB6
                           02F4
                   OFD2
                           07F4
                                           BEEP: GOTO WEITITEM
40
                   OFD9
                           02F4
                                   LEFTAR:
                   OFD9
                           02F4
                                            IF ITEN = 0 THEN GOTO NEITITEN
                           02F4
                   OFDE
                                            GOSUB ITEMBOXERASE
                   OFEE
                           02F4
                                            ITEM = ITEM - 1
                   OFF4
                           02F4
45
                                            RETIVEN DTOO
                    1003
                           02F4
                   1007
                           02F4
                                    RIGHTAR:
                           02F4
                    1007
                                            IF ITEM . MAXITEM THEN SOTO MEXITIES
                    100C
                           02F4
                    101F
                                            GOSUB ITEMBOTERASE
                           02F4
50
                                            ITEM = ITEM + 1
                    1025
                           02F4
                                            GOTO KENITER
                    1034
                           02F4
                    1038
                            02F4
                                    MENU1:
                    1038
                           02F4
                    103D
                            02F4
                                            RESTORE MIL
55
                    .1044
                            02F4
                                            RETURN
                    1048
                            02F4
```

15

20

25

RESTORE MNZ

1048

IO4D

02F4

02F4

KENUZ: